## **PRODUCT INFORMATION & MANUAL**

**Human CD25/IL-2 R $\alpha$  Valukine™ ELISA**

**VAL122**

For the quantitative determination of natural and recombinant  
human CD25/IL-2 R $\alpha$  concentrations

For research use only.  
Not for diagnostic or therapeutic procedures.

**Bio-Techne China Co. Ltd**

P: +86 (21) 52380373 P: 8009881270 F: +86 (21) 52381001  
**info.cn@bio-techne.com**

Please refer to the kit label for expiry date.  
Novus kits are guaranteed for 3 months from date of receipt

Version202209.4

## TABLE OF CONTENTS

I. BACKGROUND .....	2
II. OVERVIEW .....	3
III. ADVANTAGES .....	4
IV. EXPERIMENT .....	6
V. KIT COMPONENTS AND STORAGE .....	7
VI. PREPARATION .....	9
VII. ASSAY PROCEDURE .....	11
VIII. REFERENCES .....	13

## I. BACKGROUND

The biological activities of IL-2 are mediated by its binding to a multi-molecular cellular receptor complex. A model of the IL-2 receptor complex (1-5) would describe the high affinity receptor as an  $\alpha\beta\gamma$  trimer, in which all three chains are in contact with the ligand, it transduces IL-2 signals.

IL-2 R $\alpha$  (also known as Tac antigen and as CD25) is a 55kDa transmembrane glycoprotein composed of 351 amino acids with only 13 located on the cytoplasmic side of the membrane. IL-2 R $\alpha$  lacks structural features characteristic of members of the cytokine receptor superfamily. By itself, IL-2 R $\alpha$  binds IL-2 with low affinity (6-8).

A soluble form of IL-2 R $\alpha$  appears in serum, concomitant with its increased expression on cells (9-10). The function of the soluble IL-2 R $\alpha$  is unclear, since it would be expected to be a poor inhibitor of IL-2 because of its low binding affinity. In any case, increased levels of the soluble IL-2 R $\alpha$  in biological fluids reportedly correlate with increased T and B cell activation and immune system activation. Results of a number of studies suggest a correlation of levels of IL-2 R $\alpha$  in serum with the onset of rejection episodes in allograft recipients (9, 11-12), with activity of autoimmune diseases such as rheumatoid arthritis and systemic lupus erythematosus (SLE) (13) and with the course of some leukemias and lymphomas (14-18).

## II. OVERVIEW

### A. PRINCIPLE OF THE ASSAY

This assay employs the quantitative sandwich enzyme immunoassay technique. An antibody specific for human CD25/IL-2 R $\alpha$  has been pre-coated onto a microplate. Standards and samples are pipetted into the wells and any CD25/IL-2 R $\alpha$  present is bound by the immobilized antibody. After washing away any unbound substances, a biotin-linked detect antibody specific for human CD25/IL-2 R $\alpha$  are pipetted into the wells. After washing away any unbound substances, streptavidin-HRP is added. Following a wash to remove any unbound antibody-enzyme reagent, a substrate solution is added to the wells and color develops in proportion to the amount of CD25/IL-2 R $\alpha$  bound in the initial step. The color development is stopped, and the intensity of the color is measured.

### B. LIMITATIONS OF THE PROCEDURE

- **FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.**
- This kit is suitable for cell culture supernate and serum.
- The kit should not be used beyond the expiration date on the kit label.
- Do not mix or substitute reagents with those from other lots or sources.
- If samples generate values higher than the highest standard, dilute the samples with Calibrator Diluent and repeat the assay.
- Any variation in operator, pipetting technique, washing technique, incubation time or temperature, and kit age can cause variation in binding.

### III. ADVANTAGES

#### A. PRECISION

**Intra-assay Precision** (Precision within an assay)

Three samples were tested twenty times on one plate to assess intra-assay precision.

**Inter-assay Precision** (Precision between assays)

Three samples were tested in twenty separate assays to assess inter-assay precision.

Sample	Intra-assay Precision			Inter-assay Precision		
	1	2	3	1	2	3
Mean (pg/mL)	90.4	253.0	867.7	89.1	251.9	826.0
Standard Deviation	4.4	9.7	25.4	4.5	9.5	50.4
CV%	4.9	3.8	2.9	5.0	3.8	6.1

#### B. RECOVERY

The recovery of human CD25/IL-2 R $\alpha$  spiked to three levels throughout the range of the assay in cell culture media was evaluated. The recovery ranged from 96.3-108.6% with an average of 104.1%.

#### C. SENSITIVITY

The minimum detectable dose (MDD) of CD25/IL-2 R $\alpha$  is typically less than 9.15 pg/mL.

The MDD was determined by adding two standard deviations to the mean optical density value of twenty zero standard replicates and calculating the corresponding concentration.

#### D. CALIBRATION

This immunoassay is calibrated against a highly purified NS0-expressed recombinant human CD25/IL-2 R $\alpha$  produced at R&D Systems®.

## E. LINEARITY

To assess the linearity of the assay, different samples containing or spiked with high concentrations of CD25/IL-2 R $\alpha$  and diluted with Calibrator Diluent 1 $\times$  to produce samples with values within the dynamic range of the assay.

<b>Dilution</b>	<b>Average % of Expected</b>	<b>Range (%)</b>
1:2	100.8	89.9-107.6
1:4	94.4	83.1-105.2
1:8	91.5	82.9-104.5
1:16	88.2	84.2-103.2

## F. SAMPLE VALUES

Five serum samples were evaluated for the presence of CD25/IL-2 R $\alpha$  in this assay. All samples measured ranged from 635 to 3625.8 pg/mL with an average of 2252.2 pg/mL.

## G. SPECIFICITY

This assay recognizes both natural and recombinant human CD25/IL-2 R $\alpha$ . The following factors prepared at 50 ng/mL were assayed and exhibited no cross-reactivity or interference.


<b>Recombinant human</b>	<b>Other recombinants</b>
IL-2	mouse IL-2
IL-2 R $\gamma$	porcine IL-2
	rat IL-2

A 50 ng/mL sample of recombinant human IL-2 R $\beta$  did not exhibit any cross reactivity, but did decrease the observed reading of recombinant human CD25/IL-2 R $\alpha$  from 2.0 ng/mL to 1.6 ng/mL (20% decrease).

## IV. EXPERIMENT

### EXAMPLE STANDARD

The standard curve is provided for demonstration only. A standard curve should be generated for each set of samples assayed.
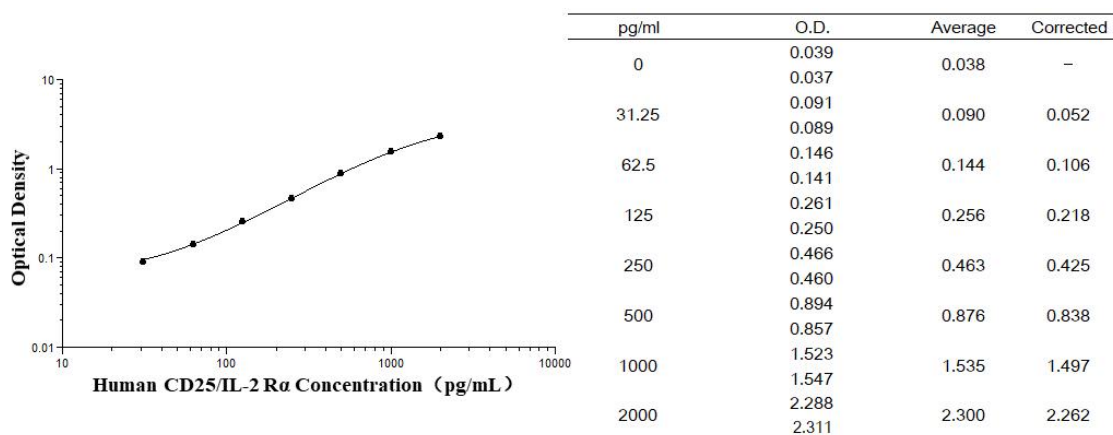
## V. KIT COMPONENTS AND STORAGE

### A. MATERIALS PROVIDED

Store the unopened kit at 2-8°C. Do not use past kit expiration date.

Parts	Description	Size
Human CD25/IL-2 R $\alpha$ Microplate	96 well polystyrene microplate (12 strips of 8 wells) coated with a mouse antibody against human CD25/IL-2 R $\alpha$ .	1 plate
Human CD25/IL-2 R $\alpha$ Standard	Recombinant human CD25/IL-2 R $\alpha$ in a buffered protein base; lyophilized. Refer to the vial label for reconstitution volume.	2 vials
Human CD25/IL-2 R $\alpha$ Detection Antibody	Biotinylated CD25/IL-2 R $\alpha$ polyclonal antibody, lyophilized. Refer to the vial label for reconstitution volume.	1 vial
Streptavidin-HRP A (200 $\times$ )	200 $\times$ Streptavidin conjugated to horseradish peroxidase.	1 vial
Reagent Diluent (10 $\times$ )	A 10 $\times$ concentrated buffered protein base used to dilute HRP.	1 vial
Calibrator Diluent (4 $\times$ )	A 4 $\times$ concentrated buffered protein base used to dilute standard and samples.	1 vial
Detection Antibody Diluent (4 $\times$ )	A buffered protein base used to dilute detection antibody.	1 vial
Wash Buffer Concentrate (25 $\times$ )	A 25 $\times$ concentrated solution of buffered surfactant with preservatives.	1 vial
TMB Substrate	TMB ELISA Substrate Solution	1 vial
Stop Solution	2 N sulfuric acid.	1 vial
Plate Sealers	Adhesive strips.	3 strips

## B. STORAGE

<b>Unopened Kit</b>	Store at 2-8 °C. Do not use past kit expiration date.	
Opened/ Reconstituted Reagents	Streptavidin-HRP A	May be stored for up to 1 month at 2-8 °C.*
	Diluted Wash Solution	
	TMB Substrate	
	Stop Solution	
	Standard	Prepare fresh for each assay.
	Detection Antibody	Aliquot and store for up to 1 month at -20 °C in a manual defrost freezer. *
	Reagent Diluent (10×)	May be stored for up to 1 month at 2-8 °C.*  Use and discard diluted Reagent Diluent (1×). Prepare fresh for each assay.
	Calibrator Diluent (4×)	May be stored for up to 1 month at 2-8 °C.*  Use and discard diluted Calibrator Diluent (1×). Prepare fresh for each assay.
	Detection Antibody Diluent (4×)	May be stored for up to 1 month at 2-8 °C.*
Microplate Wells	Return unused wells to the foil pouch containing the desiccant pack, reseal along entire edge of zip-seal. May be stored for up to 1 month at 2-8 °C.*	

\* Provided this is within the expiration date of the kit.

## C. OTHER SUPPLIES REQUIRED

- ◆ Microplate reader capable of measuring absorbance at 450 nm, with the correction wavelength set at 540 nm or 570 nm.
- ◆ Pipettes and pipette tips.
- ◆ Deionized or distilled water.
- ◆ Squirt bottle, manifold dispenser, or automated microplate washer.
- ◆ Test tubes for dilution of standards.
- ◆ 100 mL and 500 mL graduated cylinder.

## D. PRECAUTION

- ◆ Some components in this kit contain a preservative which may cause an allergic skin reaction. Avoid breathing mist.
- ◆ The Stop Solution provided with this kit is an acid solution.
- ◆ Wear protective gloves, clothing, eye, and face protection. Wash hands thoroughly after handling.

## VI. PREPARATION

### A. SAMPLE COLLECTION AND STORAGE

**Cell Culture Supernates** - Remove particulates by centrifugation and assay immediately or aliquot and store samples at  $\leq -20$  °C. Avoid repeated freeze-thaw cycles.

**Serum** - Use a serum separator tube (SST) and allow samples to clot for 30 minutes at room temperature before centrifugation for 15 minutes at 1000 x g. Remove serum and assay immediately or aliquot and store samples at  $\leq -20$  °C. Avoid repeated freeze-thaw cycles.

### B. SAMPLE PREPARATION

Cell culture supernate samples require at least a 10-fold dilution prior to the assay. A suggested 10-fold dilution is 25  $\mu$ L of sample + 225  $\mu$ L of Calibrator Diluent (1 $\times$ ).

Serum samples require at least a 4-fold dilution prior to the assay. A suggested 4-fold dilution is 50  $\mu$ L of sample + 150  $\mu$ L of Calibrator Diluent (1 $\times$ ).

### C. REAGENT PREPARATION

**Note:** *Bring all reagents to room temperature before use.*

**Wash Buffer** - If crystals have formed in the concentrate, warm to room temperature and mix gently until the crystals have completely dissolved. Dilute 20 mL of Wash Buffer Concentrate (25 $\times$ ) into deionized or distilled water to prepare 500 mL of Wash Buffer.

**Reagent Diluent (1 $\times$ )** - Add 3 mL of Reagent Diluent (10 $\times$ ) into 27 mL of deionized or distilled water to prepare 30 mL of Reagent Diluent (1 $\times$ ).

**Calibrator Diluent (1 $\times$ )** - Use deionized or distilled water to prepare Calibrator Diluent (1 $\times$ ).

**Detection Antibody Diluent (1 $\times$ )** - Use deionized or distilled water to prepare Detection Antibody Diluent (1 $\times$ ).

**Detection Antibody** - **Centrifuge briefly before opening. Reconstitution Volume refer to vial label with Reagent Diluent (1 $\times$ ).** Aliquot and store if needed. Dilute stock solution in Detection Antibody Diluent (1 $\times$ ) to the working concentration of 100 ng/mL. Prepare at least 15 minutes prior to use.


**Streptavidin-HRP A (1 $\times$ )** - **Centrifuge briefly before opening.** Dilute to the working concentration specified on the vial label using Reagent Diluent (1 $\times$ ).

**CD25/IL-2 R $\alpha$  Standard** – **Centrifuge briefly before opening. Refer to the vial label for the reconstitution volume\* using Calibrator Diluent 1 $\times$ .** This reconstitution

produces a stock solution of 4000 pg/mL. Allow the standard to sit for a minimum of 15 minutes with gentle agitation prior to making dilutions.

\*if you have any question, please seek help from our Technical Support.

**Pipette 300  $\mu$ L of Calibrator Diluent 1 $\times$  into each tube.** Use the stock solution to produce a dilution series (below). Mix each tube thoroughly before the next transfer. The 2000 pg/mL standard serves as the high standard. The **Calibrator Diluent 1 $\times$**  serves as the zero standard (0 pg/mL).
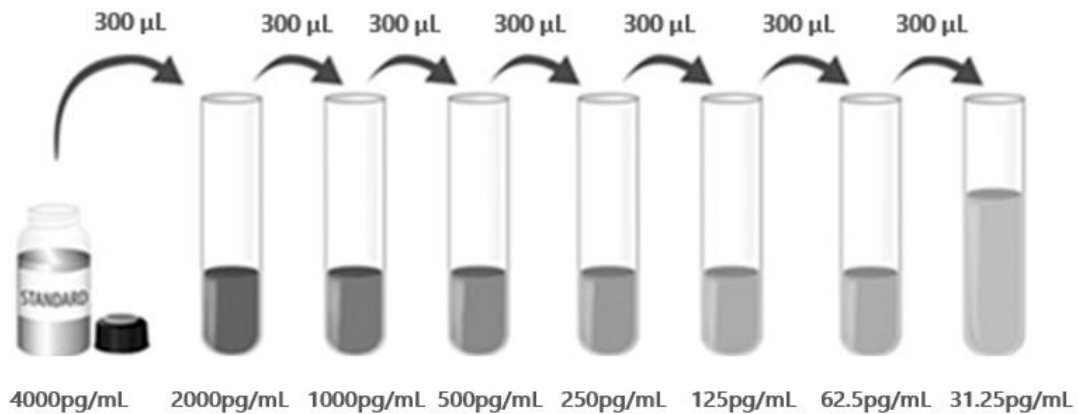
#### D. TECHNICAL HINTS

- When mixing or reconstituting protein solutions, always avoid foaming.
- To avoid cross-contamination, change pipette tips between additions of each standard level, between sample additions, and between reagent additions. Also, use separate reservoirs for each reagent.
- It is recommended that the samples be pipetted within 15 minutes.
- To ensure accurate results, proper adhesion of plate sealers during incubation steps is necessary.
- TMB Substrate should remain colorless until added to the plate. Keep Substrate Solution protected from light. Substrate Solution should change from colorless to gradations of blue.
- Stop Solution should be added to the plate in the same order as the Substrate Solution. The color developed in the wells will turn from blue to yellow upon addition of the Stop Solution. Wells that are green in color indicate that the Stop Solution has not mixed thoroughly with the Substrate Solution.

## VII. ASSAY PROCEDURE

**Note:** *Bring all reagents and samples to room temperature before use. It is recommended that all samples and standards be assayed in duplicate.*

1. Prepare all reagents and working standards as directed in the previous sections.
2. Remove excess microplate strips from the plate frame, return them to the foil pouch containing the desiccant pack, and reseal.
3. Add 100  $\mu$ L of Standard, or prepared sample per well. Cover with the adhesive strip provided. Incubate for 2 hours at room temperature. A plate layout is provided for a record of standards and samples assayed. (Samples may require dilution. See Sample Preparation section.)
4. Aspirate each well and wash, repeating the process three times for a total of four washes. Wash by filling each well with Wash Solution (400  $\mu$ L) using a squirt bottle, manifold dispenser, or autowasher. Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining Wash Buffer by aspirating or decanting. Invert the plate and blot it against clean paper towels.
5. Add 100  $\mu$ L of the Detection Antibody diluted in Detection Antibody Diluent, to each well. Cover with a new adhesive strip and incubate 2 hours at room temperature.
6. Repeat the aspiration/wash as in step 4.
7. Add 100  $\mu$ L of the working dilution of Streptavidin-HRP A to each well. Cover the plate and incubate for 20 minutes at room temperature. Avoid placing the plate in direct light.
8. Repeat the aspiration/wash as in step 4.
9. Add 100  $\mu$ L of TMB Substrate to each well. Incubate for 20 minutes at room temperature. Avoid placing the plate in direct light.
10. Add 50  $\mu$ L of Stop Solution to each well. Gently tap the plate to ensure thorough mixing.
11. Determine the optical density of each well within 10 minutes, using a microplate reader set to 450 nm. If wavelength correction is available, set to 540 nm or 570 nm. If wavelength correction is not available, subtract readings at 540 nm or 570 nm from the readings at 450 nm. This subtraction will correct for optical imperfections in the plate. Readings made directly at 450 nm without correction may be higher and less accurate.

## 12. CALCULATION OF RESULTS

Average the duplicate readings for each standard, control, and sample and subtract the average zero standard optical density (O.D.).

Create a standard curve by reducing the data using computer software capable of generating a four parameter logistic (4-PL) curve-fit. As an alternative, construct a standard curve by plotting the mean absorbance for each standard on the y-axis against the concentration on the x-axis and draw a best fit curve through the points on the graph. The data may be linearized by plotting the log of the CD25/IL-2 R $\alpha$  concentrations versus the log of the O.D. and the best fit line can be determined by regression analysis. This procedure will produce an adequate but less precise fit of the data.

If samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor.

## VIII. REFERENCES

1. Smith, K.A. (1989) *Annu. Rev. Cell Biol.* 5:397.
2. Minami, Y. et al. (1993) *Annu. Rev. Immunol.* 11:245.
3. Taniguchi, T. and Y. Minami (1993) *Cell* 73:5.
4. Voss, S.D. et al. (1994) *Blood* 83:626.
5. Leonard, W.J. et al. (1994) *Curr. Opin. Immunol.* 6:631.
6. Leonard, W.J. et al. (1984) *Nature* 311:626.
7. Nikaido, T. et al. (1984) *Nature* 311:631.
8. Cosman, D. et al. (1984) *Nature* 312:768.
9. Waldmann, T.A. (1993) *Immunol. Today* 14:264.
10. Fernandez-Botran, R. (1991) *FASEB J.* 5:2567.
11. Colvin, R.B. et al. (1987) *Clin. Immunol. Immunopathol.* 43:273.
12. Southern, J.F. et al. (1987) *J. Heart Transplant.* 5:370.
13. Wolf, R.E. et al. (1988) *Arthritis Rheum.* 31:729.
14. Pizzolo, G. et al. (1987) *Br. J. Haematol.* 67:377.
15. Wagner, D.K. et al. (1987) *J. Clin. Oncol.* 5:1262.
16. Semenzato, G. et al. (1987) *Blood* 70:396.
17. Chilosi, M. et al. (1987) *Blood* 70:1530.
18. Steis, R.G. et al. (1988) *Blood* 71:1304.

**PLATE LAYOUT**

Use this plate layout to record standards and samples assayed.

<b>1</b>											
<b>2</b>											
<b>3</b>											
<b>4</b>											
<b>5</b>											
<b>6</b>											
<b>7</b>											
<b>8</b>											
<b>9</b>											
<b>10</b>											
<b>11</b>											
<b>12</b>											
	<b>A</b>	<b>B</b>	<b>C</b>	<b>D</b>	<b>E</b>	<b>F</b>	<b>G</b>	<b>H</b>			

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