

Recombinant SARS-CoV-2 Spike S1 Subunit Fc Chimera


Catalog Number: 10622-CV

DESCRIPTION				
Source	Trichoplusia ni, T. ni (baculovirus)-derived sars-cov-2 Spike S1 Subunit protein			
	SARS-CoV-2 Spike S1 Subunit (Val16-Pro681) Accession # YP_009724390.1	IEGRMD	Human IgG ₁ (Pro100-Lys330)	
	N-terminus		C-terminus	
N-terminal Sequence Analysis	Val16			
Structure / Form	Disulfide-linked homodimer			
Predicted Molecular Mass	101 kDa			

SPECIFICATIONS		
SDS-PAGE	105-118 kDa, under reducing conditions	
Activity	Measured by its binding ability in a functional ELISA with Recombinant Human ACE-2 His-tag (Catalog # 933-ZN).	
Endotoxin Level	<0.30 EU per 1 μ g of the protein by the LAL method.	
Purity	>95%, by SDS-PAGE visualized with Silver Staining and quantitative densitometry by Coomassie® Blue Staining.	
Formulation	Lyophilized from a 0.2 µm filtered solution in PBS with Trehalose. See Certificate of Analysis for details.	

PREPARATION AND STORAGE			
Reconstitution	Reconstitute at 500 μg/mL in PBS.		
Shipping	The product is shipped at ambient temperature. Upon receipt, store it immediately at the temperature recommended below.		
Stability & Storage	Use a manual defrost freezer and avoid repeated freeze-thaw cycles.		
	 12 months from date of receipt, -20 to -70 °C as supplied. 		
	 1 month. 2 to 8 °C under sterile conditions after reconstitution. 		

- 1 month, 2 to 8 °C under sterile conditions after reconstitution.
- 3 months, -20 to -70 °C under sterile conditions after reconstitution.

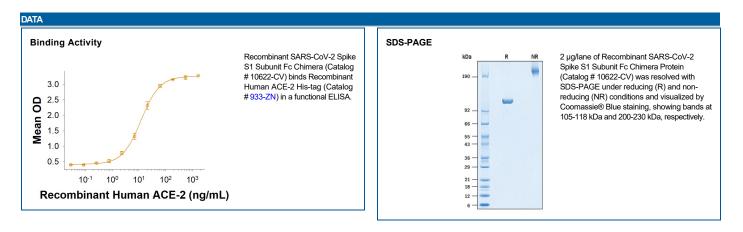

Rev. 10/13/2020 Page 1 of 2



Global bio-techne.com info@bio-techne.com techsupport@bio-techne.com TEL +1 612 379 2956 USA TEL 800 343 7475 Canada TEL 855 668 8722 China TEL +86 (21) 52380373 Europe | Middle East | Africa TEL +44 (0)1235 529449



Recombinant SARS-CoV-2 Spike S1 Subunit Fc Chimera

Catalog Number: 10622-CV

BACKGROUND

SARS-CoV-2, which causes the global pandemic coronavirus disease 2019 (Covid-19), belongs to a family of viruses known as coronaviruses that are commonly comprised of four structural proteins: Spike protein(S), Envelope protein (E), Membrane protein (M), and Nucleocapsid protein (N) (1). SARS-CoV-2 Spike Protein (S Protein) is a homotrimeric glycoprotein that mediates membrane fusion and viral entry. As with most coronaviruses, proteolytic cleavage of the SARS-CoV-2 S protein into two distinct peptides, S1 and S2 subunits, is required for activation. The S1 subunit is focused on attachment of the protein to the host receptor while the S2 subunit is involved with cell fusion (2-5). Based on structural biology studies, the receptor binding domain (RBD), located in the C-terminal region of S1, can be oriented either in the up/standing or down/lying state (6). The standing state is associated with higher pathogenicity and both SARS-CoV-1 and MERS can access this state due to the flexibility in their respective RBDs. A similar two-state structure and flexibility is found in the SARS-CoV-2 RBD (7). Based on amino acid (aa) sequence homology, the SARS-CoV-2 S1 subunit has 65% identity with SARS-CoV-1 S1 subunit, but only 22% homology with the MERS S1 subunit. The low aa sequence homology is consistent with the finding that SARS and MERS bind different cellular receptors (8). The S Protein of the SARS-CoV-2 virus, like the SARS-CoV-1 counterpart, binds Angiotensin-Converting Enzyme 2 (ACE-2), but with much higher affinity and faster binding kinetics (9). Before binding to the ACE-2 receptor, structural analysis of the S1 trimer shows that only one of the three RBD domains in the trimeric structure is in the "up" conformation. This is an unstable and transient state that passes between trimeric subunits but is nevertheless an exposed state to be targeted for neutralizing antibody therapy (10). Polyclonal antibodies to the RBD of the SARS-CoV-2 S1 subunit have been shown to inhibit interaction with the ACE-2 receptor, confirming S1 subunit especially RBD as an attractive target for vaccinations or antiviral therapy (11). There is also promising work showing that the RBD may be used to detect presence of neutralizing antibodies present in a patient's bloodstream, consistent with developed immunity after exposure to the SARS-CoV-2 virus (12). Lastly, it has been demonstrated the S Protein can invade host cells through the CD147/EMMPRIN receptor and mediate membrane fusion (13, 14).

References:

- 1. Wu, F. et al. (2020) Nature 579:265.
- 2. Tortorici, M.A. and D. Veesler (2019). Adv. Virus Res. 105:93.
- 3. Bosch, B.J. et al. (2003). J. Virol. 77:8801.
- 4. Belouzard, S. et al. (2009) Proc. Natl. Acad. Sci. 106:5871.
- 5. Millet, J.K. and G. R. Whittaker (2015) Virus Res. 202:120.
- 6. Yuan, Y. *et al.* (2017) Nat. Commun. **8**:15092.
- 7. Walls, A.C. et al. (2010) Cell 180:281.
- 8. Jiang, S. et al. (2020) Trends. Immunol. https://doi.org/10.1016/j.it.2020.03.007.
- 9. Ortega, J.T. et al. (2020) EXCLI J. 19:410.
- 10. Wrapp, D. et al. (2020) Science 367:1260.
- 11. Tai, W. et al. (2020) Cell. Mol. Immunol. https://doi.org/10.1016/j.it.2020.03.007.
- 12. Okba, N. M. A. et al. (2020). Emerg. Infect. Dis. https://doi.org/10.3201/eid2607.200841.
- 13. Wang, X. et al. (2020) https://doi.org/10.1038/s41423-020-0424-9.
- 14. Wang, K. et al. (2020) bioRxiv https://www.biorxiv.org/content/10.1101/2020.03.14.988345v1.

Rev. 10/13/2020 Page 2 of 2



Global bio-techne.com info@bio-techne.com techsupport@bio-techne.com TEL +1 612 379 2956 USA TEL 800 343 7475 Canada TEL 855 668 8722 China TEL +86 (21) 52380373 Europe | Middle East | Africa TEL +44 (0)1235 529449