

Source	Human embryonic kidney cell, HEK293-derived sars-cov-2 Spike protein			
	SARS-CoV-2 B.1.1.529 Spike (Val16 - Lys1211) (Ala67Val, His69del, Val70del, Thr95lle, Gly142Asp, Val143del, Tyr144del, Tyr145del, Asn211del, Leu212lle, ins214Glu-Pro-Glu, Gly339Asp, Ser371Leu, Ser373Pro, Ser375Phe, Lys417Asn, Asn440Lys, Gly446Ser, Ser477Asn, Thr478Lys, Glu484Ala, Gln493Arg, Gly496Ser, Gln498Arg, Asn501Tyr, Tyr505His, Thr547Lys, Asp614Gly, His655Tyr, Asn679Lys, Pro681His, Asn764Lys, Asp796Tyr, Asn856Lys, Gln954His, Asn969Lys, Leu981Phe) (Arg682Ser, Arg685Ser, Lys986Pro, Val987Pro)	6-His tag		
	Accession # YP_009724390.1 N-terminus	C-terminus		
N-terminal Sequence Analysis	Val16			
Predicted Molecular Mass	134 kDa			

SPECIFICATIONS		
SDS-PAGE	146-166 kDa, under reducing conditions.	
Activity	Measured by its binding ability in a functional ELISA with Recombinant Human ACE-2 His-tag (Catalog # 933-ZN).	
Endotoxin Level	<0.10 EU per 1 µg of the protein by the LAL method.	
Purity	>95%, by SDS-PAGE visualized with Silver Staining and quantitative densitometry by Coomassie® Blue Staining.	
Formulation	Lyophilized from a 0.2 µm filtered solution in PBS with Trehalose. See Certificate of Analysis for details.	

PREPARATION AND STORAGE		
Reconstitution	Reconstitute at 500 μg/mL in PBS.	
Shipping	The product is shipped at ambient temperature. Upon receipt, store it immediately at the temperature recommended below.	
Stability & Storage	Use a manual defrost freezer and avoid repeated freeze-thaw cycles.	
	<ul> <li>12 months from date of receipt, -20 to -70 °C as supplied.</li> </ul>	
	<ul> <li>1 month, 2 to 8 °C under sterile conditions after reconstitution.</li> </ul>	
	<ul> <li>3 months, -20 to -70 °C under sterile conditions after reconstitution.</li> </ul>	


Rev. 5/23/2024 Page 1 of 3

Bio-Techne® Global | bio-techne.com info@bio-techne.com techsupport@bio-techne.com TEL: 1.612.379.2956 USA | TEL: 800.343.7475 Canada | TEL: 855.668.8722 Europe | Middle East | Africa TEL: +44.0.1235.529449 China | info.cn@bio-techne.com TEL: 400.821.3475

## bio-techne<sup>®</sup> RDSYSTEMS


## Recombinant SARS-CoV-2 B.1.1.529 Spike His-tag

Catalog Number: 11060-CV

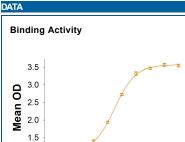


1.0 10<sup>-1</sup> 10<sup>0</sup> 10<sup>1</sup> 10<sup>2</sup> 10<sup>3</sup> Recombinant Human ACE-2 (ng/mL) Recombinant SARS-CoV-2 B.1.1.529 Spike His-tag Protein Binding Activity. Recombinant SARS-CoV-2 B.1.1.529 Spike His-tag (Catalog # 11060-CV) binds Recombinant Human ACE-2 His-tag (Catalog # 933-ZN) in a functional ELISA.



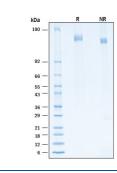

Recombinant SARS-CoV-2 B.1.1.523 Spike His-tag Protein SDS-PAGE. 2 µg/lane of Recombinant SARS-CoV-2 B.1.1.529 Spike His-tag Protein (Catalog # 11060-CV) was resolved with SDS-PAGE under reducing (R) and non-reducing (NR) conditions and visualized by Coomassie® Blue staining, showing bands at 146 - 166 kDa.

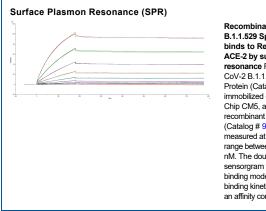

Recombinant SARS-CoV-2 B.1.1.529 Spike His-tag Protein binds to Recombinant Human ACE-2 by surface plasmon resonance Recombinant SARS-CoV-2 B.1.1.529 Spike His-tag Protein (Catalog #11060-CV) was immobilized on a Biacore Sensor Chip CM5, and binding to recombinant human ACE-2 (Catalog # 933-ZN) was measured at a concentration range between 0.021 nM and 53.1 nM. The double-referenced sensorgram was fit to a 1:1 binding model to determine the binding kinetics and affinity, with an affinity constant of Kd=3.63 nM.

Rev. 5/23/2024 Page 2 of 3

Global | bio-techne.com info@bio-techne.com techsupport@bio-techne.com TEL: 1.612.379.2956 USA | TEL: 800.343.7475 Canada | TEL: 855.668.8722 Europe | Middle East | Africa TEL: +44.0.1235.529449 China | info.cn@bio-techne.com TEL: 400.821.3475

# biotechne

## Recombinant SARS-CoV-2 B.1.1.529 Spike His-tag

**R**Dsystems

### Catalog Number: 11060-CV

#### BACKGROUND

SARS-CoV-2, which causes the global pandemic coronavirus disease 2019 (Covid-19), belongs to a family of viruses known as coronaviruses that also include MERS-CoV and SARS-CoV-1. Coronaviruses are commonly comprised of four structural proteins: Spike protein (S), Envelope protein (E), Membrane protein (M) and Nucleocapsid protein (N) (1). The SARS-CoV-2 S protein is a glycoprotein that mediates membrane fusion and viral entry. The S protein is homotrimeric, with each ~180-kDa monomer consisting of two subunits, S1 and S2 (2). In SARS-CoV-2, as with most coronaviruses, proteolytic cleavage of the S protein into S1 and S2 subunits is required for activation. The S1 subunit is focused on attachment of the protein to the host receptor while the S2 subunit is involved with cell fusion (3-5). The S protein of SARS-CoV-2 shares 75% and 29% aa sequence identity with S protein of SARS-CoV-1 and MERS, respectively. The S Protein of the SARS-CoV-2 virus, like the SARS-CoV-1 counterpart, binds a metallopeptidase, Angiotensin-Converting Enzyme 2 (ACE-2), but with much higher affinity and faster binding kinetics through the receptor binding domain (RBD) located in the C-terminal region of S1 subunit (6). It has been demonstrated that the S Protein can invade host cells through the CD147/EMMPRIN receptor and mediate membrane fusion (7, 8). Polyclonal antibodies to the RBD of the SARS-CoV-2 protein have been shown to inhibit interaction with the ACE-2 receptor, confirming RBD as an attractive target for vaccinations or antiviral therapy (9). There is also promising work showing that the RBD may be used to detect presence of neutralizing antibodies present in a patient's bloodstream, consistent with developed immunity after exposure to the SARS-CoV-2 (10). Several emerging SARS-CoV-2 genomes have been identified including the Omicron, or B.1.1.529, variant. First identified in November 2021 in South Africa, the Omicron variant quickly became the predominant SARS-CoV-2 variant and is considered a variant of concern (VOC). The Omicron variant contains 32 mutations in the S protein, 3 to 4 times more than in other SARS-CoV-2 variants, that potentially affect viral fitness and transmissibility (11). Of these mutations, 15 are located in the RBD domain and allow the Omicron variant to bind ACE-2 with greater affinity and, potentially, increased transmissibility (11, 12). Several additional mutations throughout the S protein have been shown or are predicted to enhance spike cleavage and could aid transmission (13-15). The study of the Omicron variant's impact on immune escape and reduced neutralization activity to monoclonal antibodies along with an increased risk of reinfection, even among vaccinated individuals, remains ongoing (16).

#### References:

- 1. Wu, F. et al. (2020) Nature 579:265.
- 2. Tortorici, M.A. and D. Veesler (2019) Adv. Virus Res. 105:93.
- 3. Bosch, B.J. et al. (2003) J. Virol. 77:8801.
- 4. Belouzard, S. et al. (2009) Proc. Natl. Acad. Sci. 106:5871.
- 5. Millet, J.K. and G.R. Whittaker (2015) Virus Res. 202:120.
- 6. Ortega, J.T. *et al.* (2020) EXCLI J. **19**:410.
- 7. Wang, K. et al. (2020) bioRxiv https://www.biorxiv.org/content/10.1101/2020.03.14.988345v1.
- 8. Isabel, S. *et al.* (2020) Sci Rep. 10, 14031. <u>https://doi.org/10.1038/s41598-020-70827-z</u>.
- 9. Tai, W. et al. (2020) Cell. Mol. Immunol. 17:613.
- 10. Okba, N. M. A. et al. (2020). Emerg. Infect. Dis. https://doi.org/10.3201/eid2607.200841.
- 11. Shah, M. and Woo, H.G. (2021) bioRxiv https://doi.org/10.1101/2021.12.04.471200.
- 12. Lupala, C.S. et al. (2021) bioRxiv https://doi.org/10.1101/2021.12.10.472102.
- 13. Zhang, L. et al. (2020) Nat Commun. 11:6013.
- 14. Lasek-Nesselquist, E. et al. (2021) medRxiv https://doi.org/10.1101/2021.03.10.21253285.
- 15. Scheepers, C. et al. (2021) medRxiv https://doi.org/10.1101/2021.08.20.21262342.
- 16. Callaway, E. and Ledford, H. (2021) Nature 600:197.

Rev. 5/23/2024 Page 3 of 3