

DESCRIPTION

Source	Mouse myeloma cell line, NS0-derived human HGF protein Gln ₃₂ -Ser728 Accession # P14210 Manufactured and tested under cGMP guidelines.
N-terminal Sequence Analysis	Amino acid sequencing of the alpha chain was blocked, suggesting it is consistent with Gln ₃₂ Gln ₃₂ -Arg-Lys-Arg-Arg-Asn-Thr-Ile-His-Glu (α-chain) Val ₄₉₅ -Val-Asn-Gly-Ile-Pro-Thr-Arg-Thr-Asn (β-chain)
Structure / Form	Disulfide-linked heterodimer
Predicted Molecular Mass	α chain: 53.7 kDa β chain: 26 kDa

SPECIFICATIONS


SDS-PAGE	65-75 kDa, non-reducing conditions 60-65 kDa, and 30-40 kDa, reducing conditions
Activity	Measured by its ability to induce IL-11 secretion by Saos-2 human osteosarcoma cells. Hjertner, O. <i>et al.</i> (1999) Blood 94 :3883. The ED ₅₀ for this effect is 0.1-0.5 ng/mL. The specific activity of recombinant human HGF is >1.0 x 10 ⁶ IU/mg, which is calibrated against the human HGF WHO International Standard (NIBSC code: 96/564).
Endotoxin Level	<0.10 EU per 1 µg of the protein by the LAL method.
Purity	>95%, by SDS-PAGE with silver staining, under reducing conditions.
Host Cell Protein	< 5 ng per µg of protein when tested by ELISA.
Mycoplasma	Negative when tested in a ribosomal RNA hybridization assay.
Formulation	Lyophilized from a 0.2 µm filtered solution in PBS. See Certificate of Analysis for details.

PREPARATION AND STORAGE

Reconstitution	Reconstitute at 100 µg/mL in PBS.
Shipping	The product is shipped with polar packs. Upon receipt, store it immediately at the temperature recommended below.
Stability & Storage	Use a manual defrost freezer and avoid repeated freeze-thaw cycles. <ul style="list-style-type: none"> • A minimum of 12 months when stored at ≤ -20 °C as supplied. Refer to lot specific COA for the Use by Date. • 1 month, 2 to 8 °C under sterile conditions after reconstitution. • 3 months, ≤ -20 °C under sterile conditions after reconstitution.


DATA

Bioactivity

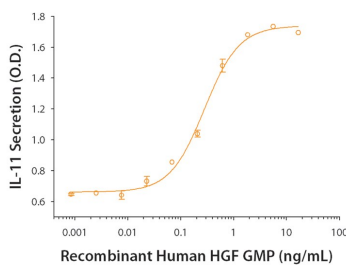

Recombinant Human HGF (NS0-expressed) GMP Protein Bioactivity Recombinant Human HGF GMP induces IL-11 secretion by Saos-2 human osteosarcoma cells. The ED₅₀ for this effect is 0.1-0.5 ng/mL.

SDS-PAGE

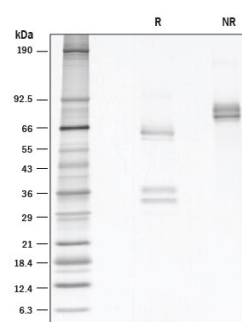

Recombinant Human HGF (NS0-expressed) GMP Protein SDS-PAGE 1 µg/lane of GMP-grade Recombinant Human HGF (Catalog # 294-GMP) was resolved with SDS-PAGE under reducing (R) and non-reducing (NR) conditions and visualized by silver staining, showing R bands at 34-37 and 60-64 kDa and NR bands at 76-83 kDa, respectively.

BACKGROUND

HGF, also known as scatter factor and hepatopoietin A, is a pleiotropic protein in the plasminogen subfamily of S1 peptidases. It is a multidomain molecule that includes an N-terminal PAN/APPLE-like domain, four Kringle domains, and a serine proteinase-like domain that has no detectable protease activity (1-5). Human HGF is secreted as an inactive 728 amino acid (aa) single chain propeptide. It is cleaved after the fourth Kringle domain by a serine protease to form bioactive disulfide-linked HGF with a 60 kDa α and 30 kDa β chain. Alternate splicing generates human HGF isoforms that lack the proteinase-like domain and different numbers of the Kringle domains. Human HGF shares 91%-94% aa sequence identity with bovine, canine, feline, mouse, and rat HGF. HGF binds heparan-sulfate proteoglycans and the widely expressed receptor tyrosine kinase, HGF R/c-MET (6, 7). HGF-dependent c-MET activation is implicated in the development of many human cancers (8). HGF regulates epithelial morphogenesis by inducing cell scattering and branching tubulogenesis (9, 10). HGF induces the up-regulation of integrin $\alpha 2 \beta 1$ in epithelial cells by a selective increase in $\alpha 2$ gene transcription (11). This integrin serves as a collagen I receptor, and its blockade disrupts epithelial cell branching tubulogenesis (11, 12). HGF can also alter epithelium morphology by the induction of nectin-1 α ectodomain shedding, an adhesion protein component of adherens junctions (13). In the thyroid, HGF induces the proliferation, motility, and loss of differentiation markers of thyrocytes and inhibits TSH-stimulated iodine uptake (14). HGF promotes the motility of cardiac stem cells in damaged myocardium (15).

References:

1. Karihaloo, A. *et al.* (2005) *Nephron Exp. Nephrol.* **100**:e40.
2. Hammond, D.E. *et al.* (2004) *Curr. Top. Microbiol. Immunol.* **286**:21.
3. Rosario, M. and W. Birchmeier (2004) *Dev. Cell* **7**:3.
4. Lesk, A.M. and W.D. Fordham (1996) *J. Mol. Biol.* **258**:501.
5. Nakamura, T. *et al.* (1989) *Nature* **342**:440.
6. Mizuno, K., *et al.* (1994) *J. Biol. Chem.* **269**:1131.
7. Gheradi, E. *et al.* (2003) *Proc. Natl. Acad. Sci.* **100**:12039.
8. Corso, S. *et al.* (2005) *Trends Mol. Med.* **11**:284.
9. Maeshima, A. *et al.* (2000) *Kid. Int.* **58**:1511.
10. Montesano, R. *et al.* (1991) *Cell* **67**:901.
11. Chiu, S-J. *et al.* (2002) *J. Biomed. Sci.* **9**:261.
12. Saelman, E.U.M. *et al.* (1995) *J. Cell Sci.* **108**:3531.
13. Tanaka, Y. *et al.* (2002) *Biochem. Biophys. Res. Commun.* **299**:472.
14. Mineo, R. *et al.* (1994) *Endocrinology* **145**:4355.
15. Urbanek, K. *et al.* (2005) *Circ. Res.* **97**:663.

MANUFACTURING SPECIFICATIONS

GMP Proteins

R&D Systems, a Bio-Techne Brand's GMP proteins are produced according to relevant sections of the following documents: WHO TRS, No. 822, 1992 Annex 1, Good Manufacturing Practices for Biological Products; USP Chapter 1043, Ancillary Materials for Cell, Gene and Tissue-Engineered Products and USP Chapter 92, Growth Factors and Cytokines Used in Cell Therapy Manufacturing.

R&D Systems' quality focus includes:

- Manufactured and tested under an ISO 9001:2015 and ISO 13485:2016 certified quality system
- Documented processes and QA control of documentation and process changes
- Personnel training programs
- Raw material testing and vendor qualification/monitoring
- Fully validated equipment, processes and test methods
- Equipment calibration schedules using a computerized calibration program
- Facility maintenance, safety programs and pest control
- Material review process for variances
- Monitoring of stability over product shelf-life

R&D Systems strives to provide our customers with the analytical characteristics of each product so that customers may determine whether our products are appropriate for their research. The Certificate of Analysis provided contains the following lot specific information:

- N-terminal amino acid analysis, SDS-PAGE analysis, and endotoxin level (as determined by LAL assay) performed on each bulk QC lot, not on individual bottlings of each QC lot
- Post-bottling lot-specific bioassay results (compliance with an established range) and results of microbial testing according to USP
- Host Cell Protein testing performed by ELISA
- Mycoplasma testing by ribosomal RNA hybridization assay

Additional testing and documentation requested by the customer can be arranged at an additional cost.

Production records and facilities are available for examination by appropriate personnel on-site at R&D Systems in Minneapolis, Minnesota USA.

R&D Systems sells GMP grade products for preclinical or clinical *ex vivo* use. They are not for *in vivo* use. Please read the following End User Terms prior to using this product.

PRODUCT SPECIFIC NOTICES

END USER TERMS OF USE OF PRODUCT

The following terms are offered to you upon your acceptance of these End User Terms of Use of Product. By using this product, you indicate your acknowledgment and agreement to these End User Terms of Use of Product. If you do not agree to be bound by and comply with all of the provisions of these End User Terms of Use of Product, you should contact your supplier of the product and make arrangements to return the product.

We suggest you print and retain a copy of these End User Terms of Use of Product for your records.

The End User is aware that R&D Systems, Inc. sells GMP products for preclinical or clinical *ex vivo* use and not for *in vivo* use. The End User further agrees, as a condition of the sale of R&D Systems' GMP products that: a) the End User will not use this GMP Product in any procedure wherein the product may be directly or indirectly administered to humans, unless the End User has obtained, or prior to their use will have obtained, an Investigational New Drug (IND) exemption from the FDA and will use the product only in accordance with the protocols of such IND and of the Institutional Review Board overseeing the proposed research, or b) the End User will use the products outside of the United States in accordance with the protocols of research approved by the Institutional Review Board or authorized ethics committee and regulatory agencies to which the End User is subject to in their territory.

R&D Systems, Inc. has the right, at its sole discretion, to modify, add or remove any terms or conditions of these End User Terms of Use without notice or liability to you. Any changes to these End User Terms of Use are effective immediately following the printing of such changes on this product insert. The most recent version of these End User Terms of Use of Product may be found at: RnDSystems.com/Legal.

You agree to review these End User Terms of Use of Product to ensure any subsequent use by you of R&D Systems' GMP Products following changes to these End User Terms of Use of Product constitutes your acceptance of all such changes.

TERMS AND CONDITIONS

The following limitation applies to R&D Systems' warranty and liability for damages: All products are warranted to meet R&D Systems' published specifications when used under normal laboratory conditions.

R&D SYSTEMS DOES NOT MAKE ANY OTHER WARRANTY OR REPRESENTATION WHATSOEVER, WHETHER EXPRESS OR IMPLIED, WITH RESPECT TO ITS PRODUCTS. IN PARTICULAR, R&D SYSTEMS DOES NOT MAKE ANY WARRANTY OF SUITABILITY, NONINFRINGEMENT, MERCHANTABILITY OR FITNESS FOR ANY PARTICULAR PURPOSE.

NOTWITHSTANDING ANY OTHER PROVISIONS OF THESE TERMS AND/OR ANY OTHER AGREEMENT BETWEEN R&D SYSTEMS AND PURCHASER FOR THE PURCHASE OF THE PRODUCTS, R&D SYSTEMS' TOTAL LIABILITY TO PURCHASER ARISING FROM OR IN RELATION TO THESE TERMS, AN AGREEMENT BETWEEN THE PARTIES OR THE PRODUCTS, WHETHER ARISING IN CONTRACT, TORT OR OTHERWISE SHALL BE LIMITED TO THE TOTAL AMOUNT PAID BY PURCHASER TO R&D SYSTEMS FOR THE APPLICABLE PRODUCTS. IN NO EVENT WILL R&D SYSTEMS BE LIABLE FOR THE COST OF PROCUREMENT OF SUBSTITUTE GOODS.

Full details of R&D Systems' Terms and Conditions of Sale can be found online at: RnDSystems.com/Legal.