


Human NCAM-1/CD56 Biotinylated Antibody

Monoclonal Rabbit IgG Clone # 2524C Catalog Number: FAB24086B 100 Tests

Species Reactivity	Human
Specificity	Detects human NCAM-1/CD56 in direct ELISAs.
Source	Monoclonal Rabbit IgG Clone # 2524C
Purification	Protein A or G purified from hybridoma culture supernatant
Immunogen	Mouse myeloma cell line, NS0-derived human NCAM -1/CD56 Leu20-Pro603, Glu636-Asn741 Accession # NP_001070150
Conjugate	Biotin Excitation Wavelength: N/A nm Emission Wavelength: N/A nm
Formulation	Supplied in a saline solution containing BSA and Sodium Azide. See Certificate of Analysis for details.
	*Contains <0.1% Sodium Azide, which is not hazardous at this concentration according to GHS classifications. Refer to the Safety Data Shee (SDS) for additional information and handling instructions.

APPLICATIONS Please Note: Optimal dilutions should be determined by each laboratory for each application. General Protocols are available in the Technical Information section on our website. Recommended Concentration Flow Cytometry 10 µL/10⁶ cells See Below

DATA

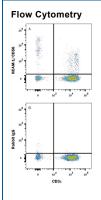

Detection of NCAM-1/CD56 in Human Peripheral Blood by Flow Cytometry. Human peripheral blood lymphocytes were stained with Mouse Anti-Human CD3 epsilon APC-conjugated Monoclonal Antibody (Catalog # FAB100A) and either (A) Rabbit Anti-Human NCAM-1/CD56 Biotinylated Monoclonal Antibody (Catalog # FAB24086B) or (B) Irrelevant antibody followed by PE-conjugated streptavidin (Catalog # F0040). Staining was performed using our Staining Membrane-associated Proteins.

PREPARATION AND STORAGE	
Shipping	The product is shipped with polar packs. Upon receipt, store it immediately at the temperature recommended below.
Stability & Storage	Do not freeze. • 12 months from date of receipt, 2 to 8 °C as supplied.





Human NCAM-1/CD56 Biotinylated Antibody

Monoclonal Rabbit IgG Clone # 2524C Catalog Number: FAB24086B

100 Tests

BACKGROUND

Neural cell adhesion molecule 1 (NCAM-1) is a multifunctional member of the Ig superfamily. It belongs to a family of membrane-bound glycoproteins that are involved in Ca⁺⁺ independent cell matrix and homophilic or heterophilic cell-cell interactions. NCAM-1 specifically binds to heparan sulfate proteoglycans (1), the extracellular matrix protein agrin (2), and several chondroitin sulfate proteoglycans that include neurocan and phosphocan (3). There are three main forms of human NCAM-1 that arise by alternate splicing. These are designated NCAM-120/NCAM-1 (761 amino acids [aa]), NCAM-140 (848 aa), and NCAM-180 (1120 aa). NCAM-120 is GPI-linked, while NCAM-140 and NCAM-180 are type I transmembrane glycoproteins (4 - 6). Additional alternate splicing adds considerable diversity to all three forms, and extracellular proteolytic processing is possible for NCAM-180 (7 - 8). NCAM-1 is synthesized as a 761 aa preproprecursor that contains a 19 aa signal sequence, a 722 aa GPI-linked mature region, and a 20 aa C-terminal prosegment (4). The molecule contains five C-2 type Ig-like domains and two fibronectin type-III domains. Human to mouse, NCAM-1 is 93% aa identical. NCAM-1 appears to be highly sialylated. The polysialyation of NCAM-1 reduces its adhesive property and increases its neurite outgrowth promoting features (9). NCAM-1 in the adult brain shows a decline of sialylation relative to earlier developmental periods. In regions that retain a high degree of neuronal plasticity, however, the adult brain continues to express polysialylation-NCAM-1, suggesting sialylation of NCAM-1 is involved in regenerative processes and synaptic plasticity (10 - 13).

References:

- 1. Burg, M.A. et al. (1995) J. Neurosci. Res. 41:49.
- 2. Storms, S.D. and U. Rutishauser (1998) J Biol. Chem. 273:27124.
- 3. Margolis, R.K. et al. (1996) Perspect. Dev. Neurobiol. 3:273.
- 4. Dickson, G. et al. (1987) Cell 50:1119.
- 5. Lanier, L.L. et al. (1991) J. Immunol. 146:4421.
- 6. Hemperly, J.J. et al. (1990) J. Mol. Neurosci. 2:71.
- 7. Rutishauser, U.and C. Goridis (1986) Trends Genet. 2:72.
- 8. Vawter, M.P. et al. (2001) Exp. Neurol. 172:29.
- 9. Rutihauser, U. (1990) Adv. Exp. Med. Biol. 265:179.
- 10. Becker, C.G. et al. (1996) J. Neurosci. Res. 45:143.
- 11. Doherty, P. et al. (1995) J. Neurobiol. 26:437.
- 12. Eckardt, M. et al. (2000) J. Neurosci. 20:5234.
- 13. Muller, D. et al. (1996) Neuron 17:413.

