biotechne

SARS-CoV-2 Spike RBD Antibody

Monoclonal Mouse IgG_{2B} Clone # 1049865 Catalog Number: MAB11290


RDSYSTEMS

DESCRIPTION		
Species Reactivity	SARS-CoV-2	
Specificity	Detects SARS-CoV2 SARS-CoV2 spike in direct ELISA.	
Source	Monoclonal Mouse IgG _{2B} Clone # 1049865	
Purification	Protein A or G purified from hybridoma culture supernatant	
Immunogen	Human embryonic kidney cell, HEK293-derived SARS-CoV2 SARS-CoV2 spike Arg319-Phe541 Accession # P0DTC2	
Formulation	Lyophilized from a 0.2 µm filtered solution in PBS with Trehalose. See Certificate of Analysis for details. *Small pack size (-SP) is supplied either lyophilized or as a 0.2 µm filtered solution in PBS.	

APPLICATIONS

Please Note: Optimal dilutions should be determined by each laboratory for each application. General Protocols are available in the Technical Information section on our website.			
Neutralization	In a functional ELISA, 0.0500 - 0.700 μg/mL of this antibody will block 50% of the binding of 50 ng/mLRecombinant Human ACE-2 (Catalog # 933-ZN) to Recombinant SARS-CoV-2 E484K N501Y Spike RBD His-tag (Catalog # 10788-CV) immobilized at 0.5 ug/mL (100 μL/well).		
Blockade of Receptor-ligand Interaction	In a functional flow cytometry test, 25 μg/mL of Mouse Anti-SARS2 RBD Antibody (Catalog # MAB11290) will block the binding of Recombinant SARS-CoV2- RBD-E484K N501Y protein (Catalog # 10788-CV) to HEK293 human embryonic kidney cell line transfected with recombinant human ACE-2.		

DATA

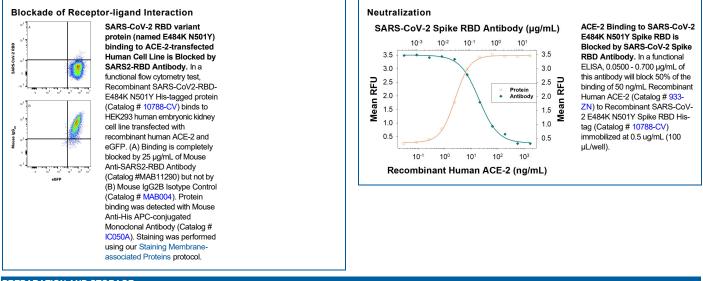

PREPARATION AND STORAGE		
Reconstitution	Reconstitute at 0.5 mg/mL in sterile PBS.	
Shipping	The product is shipped at ambient temperature. Upon receipt, store it immediately at the temperature recommended below. *Small pack size (-SP) is shipped with polar packs. Upon receipt, store it immediately at -20 to -70 °C	
Stability & Storage	Use a manual defrost freezer and avoid repeated freeze-thaw cycles.	
	 12 months from date of receipt, -20 to -70 °C as supplied. 	
	 1 month, 2 to 8 °C under sterile conditions after reconstitution. 	
	 6 months, -20 to -70 °C under sterile conditions after reconstitution. 	

Rev. 12/20/2022 Page 1 of 2



Global bio-techne.com info@bio-techne.com techsupport@bio-techne.com TEL +1 612 379 2956 USA TEL 800 343 7475 Canada TEL 855 668 8722 China TEL +86 (21) 52380373 Europe | Middle East | Africa TEL +44 (0)1235 529449

biotechne

SARS-CoV-2 Spike RBD Antibody

Monoclonal Mouse IgG_{2B} Clone # 1049865 Catalog Number: MAB11290

RDSYSTEMS

BACKGROUND

SARS-CoV-2, which caused the global pandemic coronavirus disease 2019 (Covid-19), belongs to a family of viruses known as coronaviruses that are commonly comprised of four structural proteins: Spike protein(S), Envelope protein (E), Membrane protein (M), and Nucleocapsid protein (N) (1). SARS-CoV-2 Spike Protein (S Protein) is a glycoprotein that mediates membrane fusion and viral entry. The S protein is homotrimeric, with each ~180-kDa monomer consisting of two subunits, S1 and S2 (2). In SARS-CoV-2, as with most coronaviruses, proteolytic cleavage of the S protein into two distinct peptides, S1 and S2 subunits, is required for activation. The S1 subunit is focused on attachment of the protein to the host receptor while the S2 subunit is involved with cell fusion (3-5). Based on structural biology studies, the receptor binding domain (RBD), located in the C-terminal region of S1, can be oriented either in the up/standing or down/lying state (6). The standing state is associated with higher pathogenicity and both SARS-CoV-1 and MERS can access this state due to the flexibility in their respective RBDs. A similar two-state structure and flexibility is found in the SARS-CoV-2 RBD (7). Based on amino acid (aa) sequence homology, the SARS-CoV-2 S1 subunit RBD has 73% identity with the RBD of the SARS-CoV-1 S1 RBD, but only 22% homology with the MERS S1 RBD. The low aa sequence homology is consistent with the finding that SARS and MERS bind different cellular receptors (8). The S Protein of the SARS-CoV-2 virus, like the SARS-CoV-1 counterpart, binds Angiotensin-Converting Enzyme 2 (ACE2), but with much higher affinity and faster binding kinetics (9). Before binding to the ACE2 receptor, structural analysis of the S1 trimer shows that only one of the three RBD domains in the trimeric structure is in the "up" conformation. This is an unstable and transient state that passes between trimeric subunits but is nevertheless an exposed state to be targeted for neutralizing antibody therapy (10). Polyclonal antibodies to the RBD of the SARS-CoV-2 protein have been shown to inhibit interaction with the ACE2 receptor, confirming RBD as an attractive target for vaccinations or antiviral therapy (11). There is also promising work showing that the RBD may be used to detect presence of neutralizing antibodies present in a patient's bloodstream, consistent with developed immunity after exposure to the SARS-CoV-2 virus (12). Lastly, it has been demonstrated the S Protein can invade host cells through the CD147/EMMPRIN receptor and mediate membrane fusion (13, 14).

References:

- 1. Wu, F. et al. (2020) Nature 579:265.
- 2. Tortorici, M.A. and D. Veesler (2019). Adv. Virus Res. 105:93.
- 3. Bosch, B.J. et al. (2003) J. Virol. 77:8801.
- 4. Belouzard, S. et al. (2009) Proc. Natl. Acad. Sci. 106:5871.
- 5. Millet, J.K. and G. R. Whittaker (2015) Virus Res. 202:120.
- 6. Yuan, Y. *et al.* (2017) Nat. Commun. **8**:15092.
- 7. Walls, A.C. *et al.* (2010) Cell **180**:281.
- 8. Jiang, S. et al. (2020) Trends. Immunol. https://doi.org/10.1016/j.it.2020.03.007.
- 9. Ortega, J.T. *et al*. (2020) EXCLI J. **19**:410.
- 10. Wrapp, D. et al. (2020) Science 367:1260.
- 11. Tai, W. et al. (2020) Cell. Mol. Immunol. https://doi.org/10.1016/j.it.2020.03.007.
- 12. Okba, N. M. A. et al. (2020). Emerg. Infect. Dis. https://doi.org/10.3201/eid2607.200841
- 13. Wang, X. et al. (2020) <u>https://doi.org/10.1038/s41423-020-0424-9</u>.
- 14. Wang, K. et al. (2020) bioRxiv https://www.biorxiv.org/content/10.1101/2020.03.14.988345v1.

Rev. 12/20/2022 Page 2 of 2



Global bio-techne.com info@bio-techne.com techsupport@bio-techne.com TEL +1 612 379 2956 USA TEL 800 343 7475 Canada TEL 855 668 8722 China TEL +86 (21) 52380373 Europe | Middle East | Africa TEL +44 (0)1235 529449