

DESCRIPTION

Species Reactivity	Human
Specificity	Detects human Osteoadherin/OSAD in ELISA.
Source	Monoclonal Mouse IgG _{2B} Clone # 806001
Purification	Protein A or G purified from hybridoma culture supernatant
Immunogen	HEK293 human embryonic kidney cell line transfected with human Osteoadherin/OSAD Gln21-Gln421 Accession # Q99983
Formulation	Lyophilized from a 0.2 µm filtered solution in PBS with Trehalose. See Certificate of Analysis for details. *Small pack size (-SP) is supplied either lyophilized or as a 0.2 µm filtered solution in PBS.


APPLICATIONS

Please Note: Optimal dilutions should be determined by each laboratory for each application. [General Protocols](#) are available in the Technical Information section on our website.

	Recommended Concentration	Sample
Immunocytochemistry	8-25 µg/mL	See Below

DATA

Immunocytochemistry

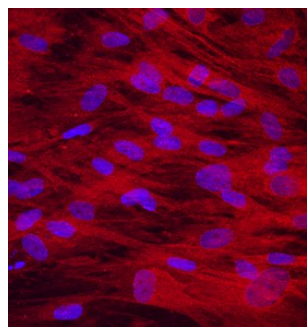

Osteoadherin/OSAD in Human Mesenchymal Stem Cells.

Osteoadherin/OSAD was detected in immersion fixed human mesenchymal stem cells differentiated into osteocytes using Mouse Anti-Human Osteoadherin/OSAD Monoclonal Antibody (Catalog # MAB2884) at 10 µg/mL for 3 hours at room temperature. Cells were stained using the NorthernLights™ 557-conjugated Anti-Mouse IgG Secondary Antibody (red; Catalog # NL007) and counterstained with DAPI (blue). Specific staining was localized to cytoplasm. View our protocol for [Fluorescent ICC Staining of Stem Cells on Coverslips](#).

PREPARATION AND STORAGE

Reconstitution	Reconstitute at 0.5 mg/mL in sterile PBS.
Shipping	The product is shipped at ambient temperature. Upon receipt, store it immediately at the temperature recommended below. *Small pack size (-SP) is shipped with polar packs. Upon receipt, store it immediately at -20 to -70 °C
Stability & Storage	Use a manual defrost freezer and avoid repeated freeze-thaw cycles. <ul style="list-style-type: none"> ● 12 months from date of receipt, -20 to -70 °C as supplied. ● 1 month, 2 to 8 °C under sterile conditions after reconstitution. ● 6 months, -20 to -70 °C under sterile conditions after reconstitution.

BACKGROUND

Osteoadherin (OSAD), also known as Osteomodulin, is an extracellular matrix keratan sulfate proteoglycan that belongs to the class II subfamily of small leucine-rich proteoglycans (SLRP). LRR motifs consist of approximately 20-30 amino acids (aa) with conserved leucine spacing, folded into a structure with one β -sheet and one α -helix (1, 2). The human OSAD cDNA encodes a 421 aa precursor that contains a 20 aa signal sequence and eleven tandem leucine rich repeats (3). Human OSAD shares 80-84% aa sequence identity with bovine, mouse, and rat OSAD. Human OSAD shares 32-35% aa sequence identity with human class II SLRPs Fibromodulin, Keratocan, Lumican, and PRELP. Bovine, mouse, and rat OSAD are expressed as 60-85 kDa molecules, although the amino acid sequence for each predicts a size of 46-47 kDa. The primary difference is due to the presence of extensive N-linked glycosylation that can also vary between tissues of the same species (4, 5). Human OSAD is expressed as an even larger 110 kDa molecule in teeth (6). OSAD contains eight sulfated tyrosine residues (4, 7) and is distinguished from other class II SLRPs by the presence of an approximately 70 aa C-terminal acidic domain (3). OSAD is expressed by fetal and adult osteoblasts but is not detectable in cartilage or tendon (3, 4, 8). In dental tissue, OSAD is expressed by odontoblasts and ameloblasts (5, 9-11) and is involved in the mineralization of bone and teeth (5, 11, 12). OSAD promotes the adhesion of osteoblasts and odontoblasts to the surrounding matrix, an interaction that is mediated by Integrin α V β 3 (4, 6).

References:

1. Matsushima, N. *et al.* (2000) *Proteins* **38**:210.
2. Kobe, B. and A.V. Kajava (2001) *Curr. Opin. Struct. Biol.* **11**:725.
3. Sommarin, Y. *et al.* (1998) *J. Biol. Chem.* **273**:16723.
4. Wendel, M. *et al.* (1998) *J. Cell Biol.* **141**:839.
5. Hultenby, P.U. *et al.* (2003) *Eur. J. Oral Sci.* **111**:128.
6. Lucchini, M. *et al.* (2004) *J. Dent. Res.* **83**:552.
7. Onnerfjord, P. *et al.* (2004) *J. Biol. Chem.* **279**:26.
8. Shen, Z. *et al.* (1999) *Matrix Biol.* **18**:533.
9. Buchaille, R. *et al.* (2000) *Bone* **27**:265.
10. Buchaille, R. *et al.* (2000) *Matrix Biol.* **19**:421.
11. Couble, M.L. *et al.* (2004) *Histochem. Cell Biol.* **121**:47.
12. Ramstad, V.E. *et al.* (2003) *Calcif. Tissue Int.* **72**:57.