

PRODUCT INFORMATION & MANUAL

Mouse IFN-β Valukine™ ELISA

VAL612

For the quantitative determination of natural and recombinant human Mouse IFN-β concentrations

For research use only.

Not for diagnostic or therapeutic procedures.

Bio-Techne China Co. Ltd
P: +86 (21) 52380373 P: 8009881270 F: +86 (21) 52381001
info.cn@bio-techne.com

Please refer to the kit label for expiry date.

Novus kits are guaranteed for 3 months from date of receipt

Version201912.1

TABLE OF CONTENTS

l.	BACKGROUND	2
II.	OVERVIEW	3
III.	ADVANTAGES	4
IV.	EXPERIMENT	6
V.	KIT COMPONENTS AND STORAGE	7
VI.	PREPARATION	. 10
VII.	ASSAY PROCEDURE	. 12
VIII	REFERENCES	14

I. BACKGROUND

Interferon beta (IFN- β), also known as fibroblast IFN, is a secreted, approximately 22 kDa member of the type I interferon family of molecules (1). Mature mouse IFN- β shares 75% and 47% amino acid sequence identity with the rat and human proteins, respectively. Fibroblasts are the major producers of IFN- β , but it can also be produced by dendritic cells, macrophages, and endothelial cells in response to pathogens (2). It is transcriptionally regulated by TRAF3, IRF3, IRF7, and NF- κ B (3, 4). IFN- β signals through the heterodimeric IFN- α / β receptor, which is constitutively associated with JAK1 and TYK2, and leads to the phosphorylation and activation of several Signal Transducer and Activator of Transcription (STAT) family members (1, 5-7).IFN- β can modulate immune responses by inducing the expression of the anti-inflammatory cytokine IL-10, promoting the generation of regulatory T cells, and inhibiting inflammasome activation (8-10).

The immune modulatory role of IFN- β appears to be important for protection against a variety of autoimmune diseases (1). IFN- β -deficient mice show increased susceptibility to experimental autoimmune encephalomyelitis (EAE), a disease model of human multiple sclerosis (MS) (11). Furthermore, IFN- β has been shown to suppress the Th17 cell response in both MS and EAE and has commonly been used as a treatment for MS (12-16). IFN- β also has protective effects when given to mice with dextran-sulphate sodium (DSS)-induced colitis (17). These immunosuppressive effects may be due in part to IFN- β -dependent induction of regulatory T cells (9, 18-20).

II. OVERVIEW

A. PRINCIPLE OF THE ASSAY

This assay employs the quantitative sandwich enzyme immunoassay technique. A monoclonal antibody specific for mouse IFN- β has been pre-coated onto a microplate. Standards and samples are pipetted into the wells and any mouse IFN- β present is bound by the immobilized antibody. After washing away any unbound substances, an enzyme-linked polyclonal antibody specific for mouse IFN- β is added to the wells. Following a wash to remove any unbound antibody-enzyme reagent, a substrate solution is added to the wells and color develops in proportion to the amount of mouse IFN- β bound in the initial step. The color development is stopped and the intensity of the color is measured.

B. LIMITATIONS OF THE PROCEDURE

- FOR RESEARCH USE ONLY NOT FOR USE IN DIAGNOSTIC PROCEDURES.
- This kit is suitable for cell culture supernate and serum.
- The kit should not be used beyond the expiration date on the kit label.
- Do not mix or substitute reagents with those from other lots or sources.
- If samples generate values higher than the highest standard, dilute the samples with Calibrator Diluent and repeat the assay.
- Any variation in operator, pipetting technique, washing technique, incubation time or temperature, and kit age can cause variation in binding.

III. ADVANTAGES

A. PRECISION

Intra-assay Precision (Precision within an assay)

Three samples were tested twenty times on one plate to assess intra-assay precision.

Inter-assay Precision (Precision between assays)

Three samples were tested in twenty separate assays to assess inter-assay precision.

	Intra-assay Precision			Inter-assay Precision		
Sample	1	2	3	1	2	3
Mean (pg/mL)	21.6	103.4	435.5	23.0	103.0	419.6
Standard Deviation	1.6	6.4	15.5	1.4	4.8	16.2
CV%	7.5	6.2	3.5	6.2	4.7	3.9

B. RECOVERY

The recovery of mouse IFN- β spiked to different levels throughout the range of the assay in cell culture media was evaluated. The recovery ranged from 102.9-116.6% with an average of 109.0%.

The recovery of mouse IFN- β spiked to different levels throughout the range of the assay in mouse serum was evaluated. The recovery ranged from 103.2-111.6% with an average of 107.4%.

C. SENSITIVITY

The minimum detectable dose (MDD) of mouse IFN- β is typically less than 5.85 pg/mL.

The MDD was determined by adding two standard deviations to the mean optical density value of twenty zero standard replicates and calculating the corresponding concentration.

D. CALIBRATION

This immunoassay is calibrated against highly purified HEK293-expressed recombinant mouse IFN-β produced at R&D Systems.

E. LINEARITY

To assess the linearity of the assay, different samples were containing or spiked with high concentrations of mouse IFN- β and diluted with Calibrator Diluent (1×) to produce samples with values within the dynamic range of the assay.

Dilution	Average % of Expected	Range (%)
1:2	96.9	91.1-103.2
1:4	92.4	91.4-99.4
1:8	84.0	76.0-87.9
1:16	79.0	70.2-87.9

F. SAMPLE VALUES

Serum - Five mouse serum samples were evaluated for the presence of mouse IFN- β in this assay. All samples measured less than the lowest mouse IFN- β standard, 12.5 pg/mL.

G. SPECIFICITY

The following factors prepared at 50 ng/mL were assayed and exhibited no cross-reactivity or interference.


Recombinant mouse	Recombinant human
IFN- α / β R2/Fc Chimera	IFN-β
ΙΕΝ- γ	
IFN- κ	
IL-6	
Limitin	

Recombinant mouse IFN- α/β R1 does not cross react in the assay does interfere at concentrations> 625 pg/mL.

IV. EXPERIMENT

EXAMPLE STANDARD

The standard curve is provided for demonstration only. A standard curve should be generated for each set of samples assayed.

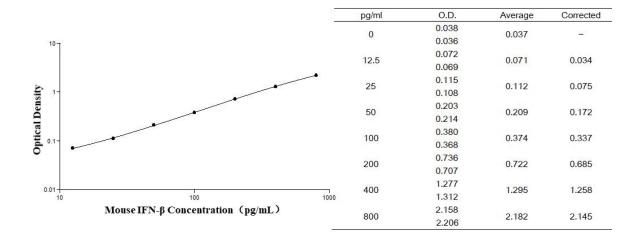

V. KIT COMPONENTS AND STORAGE

A. MATERIALS PROVIDED

Parts	Description	
Mouse IFN-β Microplate	96 well polystyrene microplate (12 strips of 8 wells) coated with a rat antibody against mouse IFN-β.	1 plate
Mouse IFN-β Standard	Recombinant mouse IFN-β in a buffered protein base; lyophilized. Refer to the vial label for reconstitution volume.	2 vials
Mouse IFN-β Detection Antibody	Biotinylated IFN-β polyclonal antibody, lyophilized. Refer to the vial label for reconstitution volume.	1 vial
Calibrator Diluent (4×)	Concentrated buffered diluent used to dilute standard and samples.	1 vial
Streptavidin-HRP B (40×)	40× Streptavidin conjugated to horseradish peroxidase.	1 vial
Reagent Diluent (10×)	A 10× concentrated buffered protein base used to dilute Detection Antibody and HRP.	1 vial
Wash Buffer Concentrate (25×)	A 25× concentrated solution of buffered surfactant with preservatives.	1 vial
Color Reagent A	Stabilized hydrogen peroxide.	1 vial
Color Reagent B	Stabilized chromogen (tetramethylbenzidine).	1 vial
Stop Solution	2 N sulfuric acid.	1 vial
Plate Sealers	Adhesive strip.	3 strips

B. STORAGE

Unopened Kit	Store at 2-8 °C. Do not use past kit expiration date.		
	Streptavidin-HRP B		
	Diluted Wash Solution		
	Unmixed Color Reagent A	May be stored for up to 1 month at 2-8 °C.*	
	Unmixed Color Reagent B		
	Stop Solution		
	Standard	Prepare fresh for each assay.	
Opened/	Detection antibody	Aliquot and store for up to 1 month at -20 °C in a manual defrost freezer. *	
Opened/ Reconstituted Reagents	Reagent Diluent (10×)	May be stored for up to 1 month at 2-8 °C.*	
		Use and discard diluted Reagent Diluent (1×). Prepare fresh for each assay.	
		May be stored for up to 1 month at 2-8 °C.*	
	Calibrator Diluent (4×)	Use and discard diluted Calibrator Diluent (1×). Prepare fresh for each assay.	
	Microplate Wells	Return unused wells to the foil pouch containing the desiccant pack, reseal along entire edge of zip-seal. May be stored for up to 1 month at 2-8 °C.*	

^{*} Provided this is within the expiration date of the kit.

C. OTHER SUPPLIES REQUIRED

- Microplate reader capable of measuring absorbance at 450 nm, with the correction wavelength set at 540 nm or 570 nm.
- Pipettes and pipette tips.
- Deionized or distilled water.
- Squirt bottle, manifold dispenser, or automated microplate washer.
- 500 mL graduated cylinder.

D. PRECAUTION

- Some components in this kit contain a preservative which may cause an allergic skin reaction. Avoid breathing mist.
- Color Reagent B may cause skin, eye, and respiratory irritation. Avoid breathing fumes.
- The Stop Solution provided with this kit is an acid solution. Wear protective gloves, clothing, eye, and face protection. Wash hands thoroughly after handling.

VI. PREPARATION

A. SAMPLE COLLECTION AND STORAGE

Cell Culture Supernates - Remove particulates by centrifugation and assay immediately or aliquot and store samples at ≤ -20°C. Avoid repeated freeze-thaw cycles. Samples may require dilution with Calibrator Diluent (1×).

Serum - Use a serum separator tube (SST) and allow samples to clot for 30 minutes at room temperature before centrifugation for 15 minutes at 1000 x g. Remove serum and assay immediately or aliquot and store samples at \leq -20°C. Avoid repeated freeze-thaw cycles.

B. SAMPLE PREPARATION

Serum samples require a 2 fold dilution. A suggested 2-fold dilution is 100 μ L of sample + 100 μ L of Calibrator Diluent (1×).

C. REAGENT PREPARATION

Note: Bring all reagents to room temperature before use.

Wash Buffer - If crystals have formed in the concentrate, warm to room temperature and mix gently until the crystals have completely dissolved. Dilute 20 mL of Wash Buffer Concentrate (25×) into deionized or distilled water to prepare 500 mL of Wash Buffer.

Calibrator Diluent (1×) - Use deionized or distilled water to prepare Calibrator Diluent (1×).

Reagent Diluent (1x) - Use deionized or distilled water to prepare Reagent Diluent (1x).

Detection Antibody- Reconstitution Volume refer to vial label with Reagent Diluent (1×). Aliquot and store if needed. Dilute stock solution in Reagent Diluent (1×) to the working concentration of 100 ng/mL. Prepare at least 15 minutes prior to use.

Streptavidin-HRP B (1×) - Dilute to the working concentration specified on the vial label using Reagent Diluent (1×).

Substrate Solution - Color Reagent A and B should be mixed together in equal volumes within 15 minutes of use. Protect from light. 100 µL of the resultant mixture is


required per well.

Mouse IFN-β Standard – Refer to the vial label for the reconstitution volume*.

This reconstitution produces a stock solution of 8000 pg/mL. Allow the standard to sit for a minimum of 15 minutes with gentle agitation prior to making dilutions.

*if you have any question, please seek help from our Technical Support.

Pipette 900 μL of Calibrator Diluent (1×) into the 800 pg/mL tube. Pipette 500 μL of the appropriate Calibrator Diluent (1×) into the remaining tubes. Use the stock solution to produce a dilution series (below). Mix each tube thoroughly before the next transfer. The 800 pg/mL standard serves as the high standard. The Calibrator Diluent (1×) serves as the zero standard (0 pg/mL).

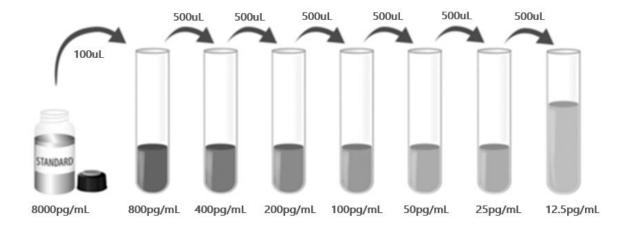

D. TECHNICAL HINTS

- When mixing or reconstituting protein solutions, always avoid foaming.
- To avoid cross-contamination, change pipette tips between additions of each standard level, between sample additions, and between reagent additions. Also, use separate reservoirs for each reagent.
- It is recommended that the samples be pipetted within 15 minutes.
- To ensure accurate results, proper adhesion of plate sealers during incubation steps is necessary.
- Substrate Solution should remain colorless until added to the plate. Keep Substrate Solution protected from light. Substrate Solution should change from colorless to gradations of blue.
- Stop Solution should be added to the plate in the same order as the Substrate Solution. The color developed in the wells will turn from blue to yellow upon addition of the Stop Solution. Wells that are green in color indicate that the Stop Solution has not mixed thoroughly with the Substrate Solution.

VII. ASSAY PROCEDURE

Note: Bring all reagents and samples to room temperature before use. It is recommended that all samples and standards be assayed in duplicate.

- 1. Prepare all reagents and working standards as directed in the previous sections.
- 2. Remove excess microplate strips from the plate frame, return them to the foil pouch containing the desiccant pack, and reseal.
- 3. Add 100 μL of Standard, or prepared sample per well. Cover with the adhesive strip provided. Incubate for 2 hours at room temperature. A plate layout is provided for a record of standards and samples assayed. (Samples may require dilution. See Sample Preparation section.)
- 4. Aspirate each well and wash, repeating the process three times for a total of four washes. Wash by filling each well with Wash Solution (400 μL) using a squirt bottle, manifold dispenser, or autowasher. Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining Wash Buffer by aspirating or decanting. Invert the plate and blot it against clean paper towels.
- 5. Add 100 μL of the Detection Antibody diluted in Reagent Diluent, to each well. Cover with a new adhesive strip and incubate 2 hours at room temperature.
- 6. Repeat the aspiration/wash as in step 4.
- 7. Add 100 μ L of the working dilution of Streptavidin-HRP B to each well. Cover the plate and incubate for 20 minutes at room temperature. Avoid placing the plate in direct light.
- 8. Repeat the aspiration/wash as in step 4.
- 9. Add 100 μL of Substrate Solution to each well. Incubate for 20 minutes at room temperature. Avoid placing the plate in direct light.
- 10. Add 50 μL of Stop Solution to each well. Gently tap the plate to ensure thorough mixing.
- 11. Determine the optical density of each well immediately, using a microplate reader set to 450 nm. If wavelength correction is available, set to 540 nm or 570 nm. If wavelength correction is not available, subtract readings at 540 nm or 570 nm from the readings at 450 nm. This subtraction will correct for optical imperfections

in the plate. Readings made directly at 450 nm without correction may be higher and less accurate.

12. CALCULATION OF RESULTS


Average the duplicate readings for each standard, control, and sample and subtract the average zero standard optical density (O.D.). Create a standard curve by reducing the data using computer software capable of generating a four parameter logistic (4-PL) curve-fit. As an alternative, construct a standard curve by plotting the mean absorbance for each standard on the y-axis against the concentration on the x-axis and draw a best fit curve through the points on the graph. The data may be linearized by plotting the log of the IFN- β concentrations versus the log of the O.D. and the best fit line can be determined by regression analysis. This procedure will produce an adequate but less precise fit of the data. If samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor.

VIII. REFERENCES

- 1. González-Navajas, J.M. et al. (2012) Nat. Rev. Immunol. 12:125.
- 2. Reder, A.T. and X. Feng (2013) Front. Immunol. 4:281.
- 3. Schafer, S.L. et al. (1998) J. Biol. Chem. 273:2714.
- 4. Häcker, H. et al. (2006) Nature 439:204.
- 5. Platanias, L.C. (2005) Nat. Rev. Immunol. 5:375.
- 6. Fasler-Kan, E. et al. (1998) Eur. J. Biochem. 254:514.
- 7. Matikainen, S. et al. (1999) Blood 93:1980.
- 8. Wang, H. et al. (2011) J. Immunol 186:675.
- 9. Qiu, H. et al. (2008) J. Immunol. 181:2092.
- 10. Guarda, G. et al. (2011) Immunity 34:213.
- 11. Teige, I. et al. (2003) J. Immunol. 170:4776.
- 12. Shinohara, M.L. et al. (2008) Immunity 29:68.
- 13. Guo, B. et al. (2008) J. Clin. Invest. 118:1680.
- Ramgolam, V.S. and S. Markovic-Plese (2010) Endocr. Metab. Immune Disord.
 Drug Targets 10:161.
- 15. Martín-Saavedra, F.M. et al. (2008) Mol. Immunol. 45:4008.
- 16. Inoue, M. and M.L. Shinohara (2013) Immunology 139:11.
- 17. Katakura, K. et al. (2005) J. Clin. Invest. 115:695.
- 18. Karaghiosoff, M. et al. (2003) Nat. Immunol. 4:471.
- 19. Vandenbark, A.A. et al. (2009) J. Neuroimmunol. 215:125.
- 20. Wang, D. et al. (2016) J. Immunol. 197:2992.

PLATE LAYOUT

Use this plate layout to record standards and samples assayed.

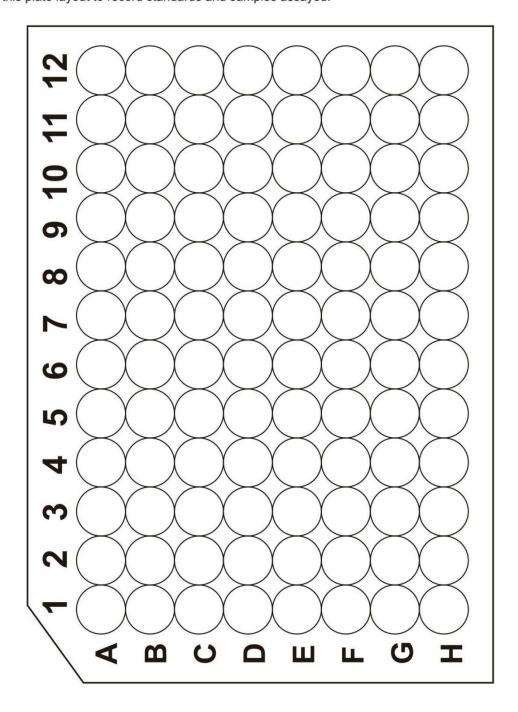



产品信息及操作手册

小鼠 IFN-β Valukine™ ELISA 试剂盒

目录号: VAL612

适用于定量检测天然和重组小鼠 IFN-β 的浓度

科研专用,不可用于临床诊断

Bio-Techne China Co. Ltd

P: +86 (21) 52380373 P: 8009881270 F: +86 (21) 52381001 info.cn@bio-techne.com

有效期详见试剂盒包装标签 Novus 试剂盒确保在你收货日期 3 个月内有效

目录

I.	背景	18
II.	概述	. 19
III.	优势	. 20
IV.	实验	. 22
V.	试剂盒组成及储存	. 23
VI.	实验前准备	. 25
VII.	操作步骤	. 27
VIII.	参考文献	28

I. 背景

干扰素β(IFN-β),也被称为成纤维细胞干扰素,是一种分泌表达的,约22 kDa的I型干扰素分子家族成员(1)。成熟的小鼠IFN-β分别与大鼠和人的蛋白有75%和47%的氨基酸序列同源性。IFN-β主要由成纤维细胞产生,但也可以由树突细胞、巨噬细胞和内皮细胞应答病原体的反应中产生(2)。它可由TRAF3、IRF3、IRF7和NF-κB转录调控(3, 4)。IFN-β通过与JAK1和TYK2有结构联系的异二聚体IFN-α/β受体发出信号,导致几个信号转导子和转录活化子(STAT)家庭成员磷酸化和激活 (1, 5-7)。IFN-β可以通过诱导抗炎细胞因子IL-10的表达,促进调节性T细胞的生成,抑制炎症小体的活化来调节免疫应答(8-10)。

IFN-β的免疫调节作用对预防多种自身免疫性疾病看来具有重要意义(1)。实验性自身免疫性脑脊髓炎(EAE)是人类多发性硬化(MS)的一种疾病模型,IFN-β缺陷小鼠显示了 EAE 易感性增加(11)。此外,IFN-β在MS和EAE 中均能抑制Th17细胞的反应,并常被用作 MS的治疗方法(12-16)。IFN-β对葡聚糖硫酸钠(DSS)诱导的小鼠结肠炎也有保护作用(17)。 这些免疫抑制作用可能部分归因于IFN-β依赖性诱导调节性T细胞(9, 18-20)。

Ⅱ. 概述

A. 检测原理

本实验采用双抗体夹心ELISA法。抗小鼠IFN-β单抗包被于微孔板上,样品和标准品中的小鼠IFN-β会与固定在板上的抗体结合,游离的成分被洗去;加入辣根过氧化酶标记的抗小鼠IFN-β多抗,未结合的抗体被洗去;溶液颜色与结合的目标蛋白成正比;加入终止液;用酶标仪测定吸光度。

B. 检测局限

- ◆ 仅供科研使用,不可用于体外诊断;
- ◆ 该试剂盒适用于细胞培养上清样本和小鼠血清样本;
- ◆ 请在试剂盒有效期内使用;
- ◆ 不同试剂盒及不同批号试剂盒的组分不能混用;
- ◆ 样本值若大于标准曲线的最高值,应将样本用标准品稀释剂(1x)稀释后重新检测;
- ◆ 检测结果的不同可由多种因素引起,包括实验人员的操作、移液器的使用方式、洗板 技术、反应时间或温度、试剂盒的储存等。

Ⅲ. 优势

A. 精确度

板内精确度(同一板内不同孔间的精确度)

已知浓度的三个样本,在同一板内分别检测20次,以确定板内精确度。

板间精确度(不同板之间的精确度)

已知浓度的三个样本,在不同板中分别检测20次,以确定板间精确度。

	板内精确度		板间精确度			
样本	1	2	3	1	2	3
平均值 (pg/mL)	21.6	103.4	435.5	23.0	103.0	419.6
标准差	1.6	6.4	15.5	1.4	4.8	16.2
CV%	7.5	6.2	3.5	6.2	4.7	3.9

B. 回收率

在细胞培养基样本中掺入检测范围内不同水平的小鼠IFN-β,测定其回收率。回收率范围在102.9-116.6%,平均回收率在109.0%。

在小鼠血清样本中掺入检测范围内不同水平的小鼠IFN-β,测定其回收率。回收率范围在103.2-111.6%,平均回收率在107.4%。

C. 灵敏度

小鼠IFN-β的最低可测剂量(MDD)一般小于5.85pg/mL。

MDD是根据20个重复的零标准品孔的吸光度值的平均值加两倍标准差计算得到的相对应浓度。

D. 校正

此ELISA试剂盒经由R&D Systems生产的HEK293表达的高纯度重组小鼠IFN-β蛋白所校正。

E. 线性

不同的样本中含有或掺入高浓度的小鼠IFN-β,然后用标准品稀释剂(1×)将样本稀释 到检测范围内,测定其线性。

稀释倍数	平均值/期待值 (%)	范围 (%)
1:2	96.9	91.1-103.2
1:4	92.4	91.4-99.4
1:8	84.0	76.0-87.9
1:16	79.0	70.2-87.9

F. 样本预值

血清样本 - 使用本试剂盒检测了5份小鼠血清样本中IFN-β的水平。5份样本的检测值均低于IFN-β最低检测限,12.5pg/mL。

G. 特异性

将以下因子配置成 50 ng/mL的浓度来检测没有观察到明显的交叉反应或干扰。


Recombinant mouse	Recombinant human
IFN-α/β R2/Fc Chimera	IFN-β
IFN-γ	
IFN-ĸ	
IL-6	
Limitin	

重组小鼠IFN-α/βR1在这个测试中没有交叉反应。但当蛋白浓度大于625 pg/mL时,可以观察到干扰。

IV. 实验

标准曲线实例

该标准曲线数据仅供参考,每次实验应绘制其对应的标准曲线。

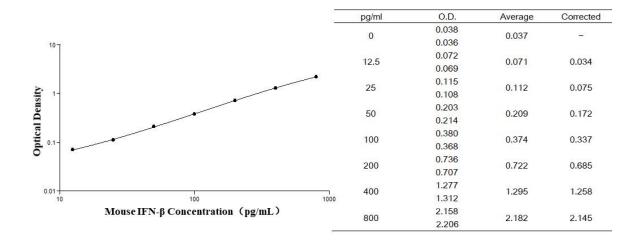

V. 试剂盒组成及储存

A. 试剂盒组成

组成	描述	规格
Mouse IFN-β Microplate	包被兔抗小鼠 IFN-β抗体的 96 孔聚苯乙烯 板,8 孔×12 条	1 块板
Mouse IFN-β Standard	标准品(冻干粉),参考瓶标签进行重溶	2 瓶
Mouse IFN-β Detection antibody	生物素化的 IFN-β检测抗体,冻干粉,参考 瓶标签进行重溶	1 瓶
Streptavidin-HRP B (40×)	40×浓缩的链霉亲和素标记的 HRP	1 瓶
Reagent Diluent (10×)	浓缩的试剂稀释液(10×)	1 瓶
Calibrator Diluent (4×)	浓缩的标准品稀释剂(4×)	1 瓶
Wash Buffer Concentrate (25×)	浓缩洗涤缓冲液(25×)	1 瓶
Color Reagent A	显色液 A	1 瓶
Color Reagent B	显色液 B	1 瓶
Stop Solution	终止液	1 瓶
Plate Sealers	封板膜	3 张

B. 试剂盒储存

未开封试剂盒	2-8℃储存;请在试剂盒有效期内使用		
	链霉亲和素-HRP B		
	洗涤缓冲液(1×)		
	显色剂 A	2-8℃储存,最多 30 天*	
	显色剂 B		
	终止液		
已打开,稀释	标准品	使用时新鲜配制*	
或重溶的试剂	检测抗体	分装, -20℃储存,最多 30 天*	
	标准品稀释剂(4×)	2-8℃储存,最多 30 天*	
	737年由中省中713 (1007)	请每次使用新鲜配制的 1×标准品稀释剂	
	试剂稀释液(10×)	2-8℃储存,最多 30 天*	
	PONTING CO	请每次使用新鲜配制的 1×试剂稀释液	
	包被的微孔板条	将未用的板条放回带有干燥剂的铝箔袋内, 密封: 2-8℃储存,最多 30 天*	

*必须在试剂盒有效期内

C. 实验所需自备试验器材

- ◆ 酶标仪(可测量450nm检测波长的吸收值及540nm或570nm校正波长的吸收值)
- ◆ 高精度加液器及一次性吸头
- ◆ 蒸馏水或去离子水
- ◆ 洗瓶 (喷瓶)、多通道洗板器或自动洗板机
- ◆ 500mL量筒

D. 注意事项

- ◆ 试剂盒中的一些组分含有防腐剂,可能引起皮肤过敏反应,避免吸入。
- ◆ Color Reagent B可能引起皮肤、眼睛和呼吸道刺激,避免吸入。
- ◆ 试剂盒中的终止液是酸性溶液,使用时请做好眼镜、手、面部及衣服的防护。使用 后请彻底洗手。

VI. 实验前准备

A. 样品收集及储存

细胞培养上清液: 颗粒物应离心去除; 立刻检测样本。样本收集后若不及时检测,需按一次使用量分装,冻存于-20℃冰箱内,避免反复冻融。样本可能需要用标准品稀释剂(1×)稀释。

血清样本:用血清分离管(SST)分离血清。使血样室温凝集30分钟,然后1000xg离心15分钟。吸取血清样本之后即刻用于检测,或者分装,-20℃贮存备用。避免反复冻融。

B. 样本准备工作

血清样本需要用标准品稀释剂 (1×) 2倍稀释后进行检测,例如: 100μL血清+100μL标准品稀释剂 (1×)。

C. 检测前准备工作

使用前请将所有试剂放置于室温

洗涤液: 从冰箱中取出的浓缩洗涤液可能有结晶,属于正常现象;放置室温,轻摇混匀,待结晶完全溶解后再配制洗涤液。可将20mL浓缩洗涤液用蒸馏水或去离子水稀释配置成500mL工作浓度的洗涤液。未用完的放回4℃。

标准品稀释剂(1×): 使用蒸馏水或去离子水稀释配置成标准品稀释剂(1×)。

试剂稀释液 (1x): 使用蒸馏水或去离子水稀释配置成试剂稀释液 (1x)。

检测抗体: 参考检测抗体瓶标签指示,用试剂稀释液(1×)将冻干粉进行重溶。再用试剂稀释液(1×)稀释至工作浓度100ng/mL,至少在使用前15分钟准备。


链霉亲和素- HRP B: 用试剂稀释液(1×)将链霉亲和素- HRP B(40×)稀释至工作浓度链霉亲和素- HRP B(1×)。

显色剂:按试验所需用量(100μL/孔)将显色液A和显色液B等体积混合,避光保存,现用现配,须在15分钟内使用。

IFN-β标准品: 冻干标准品的重溶体积请参考瓶身标签,得到浓度为8000pg/mL标准品母液。轻轻震摇至少15分钟,使其充分溶解。

*如有疑问,请咨询我们的技术支持。

800pg/mL稀释管中加入900uL标准品稀释剂(1×),其余各稀释管中加入500µL标准品稀释剂(1×)。将标准品母液参照下图做系列稀释,每管须充分混匀后再移液到下一管。800pg/mL管作标准曲线最高点,标准品稀释剂(1×)可用作标准品零点(0pg/mL)。

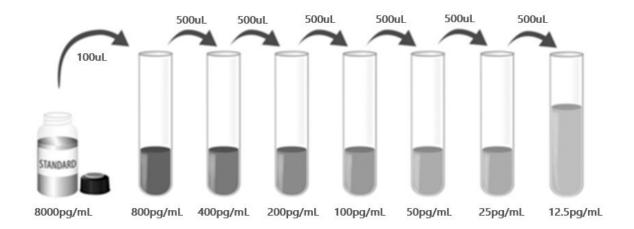

D. 技术小提示

- ◆ 当混合或重溶蛋白液时,尽量避免起沫;
- ◆ 为了避免交叉污染,配置不同浓度标准品、上样、加不同试剂都需要更换枪头。另 外不同试剂请分别使用不同的移液槽;
- ◆ 建议15分钟内完成一块板的上样;
- ◆ 每次孵育时,正确使用封板胶纸可保证结果的准确性;
- ◆ 混合后的显色底物在上板前应为无色,请避光保存;加入微孔板后,将由无色变成不同深度的蓝色;
- ◆ 终止液上板顺序应同显色底物上板顺序一致;加入终止液后,孔内颜色由蓝变黄; 若孔内有绿色,则表明孔内液体未混匀,请充分混合。

VII.操作步骤

使用前请将所有试剂和样本放置于室温,建议所有的实验样本和标准品做复孔检测

- 1. 按照上一节的说明,准备好所有需要的试剂和标准品:
- 2. 从已平衡至室温的密封袋中取出微孔板,未用的板条请放回铝箔袋内,重新封口;
- 3. 分别将不同浓度标准品,实验样本或者质控品加入相应孔中,每孔100μL。用封板膜封住反应孔,室温孵育2小时。说明书提供了一张96孔模板图,可用于记录标准品和试验样本的板内位置; (血清样本需要稀释,详情参见样本制备部分。)
- 4. 将板内液体吸去,使用洗瓶、多通道洗板器或自动洗板机洗板。每孔加洗涤液400μL,然后将板内洗涤液吸去。重复操作3次,共洗4次。每次洗板尽量吸去残留液体会有助于得到好的实验结果。最后一次洗板结束,请将板内所有液体吸干或将板倒置,在吸水纸拍干所有残留液体;
- 5. 在每个微孔内加入100μL配置好的检测抗体。用封板胶纸封住反应孔,室温孵育2 小时:
- 6. 重复第4步洗板操作:
- 7. 在每个微孔内加入100μL稀释好的链霉亲和素- HRP B工作液。用封板膜封住反应 孔,室温孵育20分钟,**注意避光**;
- 8. 重复第4步洗板操作;
- 9. 在每个微孔内加入100µL显色剂,室温孵育20分钟,注意避光;
- 10. 在每个微孔内加入50µL终止液,请轻拍微孔板,使溶液混合均匀。孔内溶液颜色会从蓝色变为黄色;
- 11. 加入终止液后30分钟内,使用酶标仪测量450nm的吸光度值,设定540nm或570nm 作为校正波长。如果没有使用双波长校正,结果准确度可能会受影响;
- 12. 计算结果:将每个标准品和样品的校正吸光度值(OD450-OD540/OD570),复孔读数取平均值,然后减去平均零标准品OD值。使用计算机软件作四参数逻辑(4-PL)曲线拟合创建标准曲线。另一种方法是,可以通过绘制标准品浓度做对数与相应OD值对数生成曲线,并通过回归分析确定最佳拟合线。通过样本的OD值,可从标准曲线上得到样本中小鼠 IFN-β的浓度。如果样品被稀释,从标准曲线读取的浓度必须乘以稀释倍数。

VIII. 参考文献

- 1. González-Navajas, J.M. et al. (2012) Nat. Rev. Immunol. 12:125.
- 2. Reder, A.T. and X. Feng (2013) Front. Immunol. 4:281.
- 3. Schafer, S.L. et al. (1998) J. Biol. Chem. 273:2714.
- 4. Häcker, H. et al. (2006) Nature 439:204.
- 5. Platanias, L.C. (2005) Nat. Rev. Immunol. 5:375.
- 6. Fasler-Kan, E. et al. (1998) Eur. J. Biochem. 254:514.
- 7. Matikainen, S. et al. (1999) Blood 93:1980.
- 8. Wang, H. et al. (2011) J. Immunol 186:675.
- 9. Qiu, H. et al. (2008) J. Immunol. 181:2092.
- 10. Guarda, G. et al. (2011) Immunity 34:213.
- 11. Teige, I. et al. (2003) J. Immunol. 170:4776.
- 12. Shinohara, M.L. et al. (2008) Immunity 29:68.
- 13. Guo, B. et al. (2008) J. Clin. Invest. 118:1680.
- Ramgolam, V.S. and S. Markovic-Plese (2010) Endocr. Metab. Immune Disord.
 Drug Targets 10:161.
- 15. Martín-Saavedra, F.M. et al. (2008) Mol. Immunol. 45:4008.
- 16. Inoue, M. and M.L. Shinohara (2013) Immunology 139:11.
- 17. Katakura, K. et al. (2005) J. Clin. Invest. 115:695.
- 18. Karaghiosoff, M. et al. (2003) Nat. Immunol. 4:471.
- 19. Vandenbark, A.A. et al. (2009) J. Neuroimmunol. 215:125.
- 20. Wang, D. et al. (2016) J. Immunol. 197:2992.

96 孔模板图

请使用 96 孔模板图来记录标准品及样本在板内的位置

