

# Product Datasheet

## Lightning-Link (R) PE-Cy5.5 Antibody Labeling Kit 761-0015

Unit Size: 600 ug

Store at -20C.

[www.novusbio.com](http://www.novusbio.com)



[technical@novusbio.com](mailto:technical@novusbio.com)

**Publications: 8**

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:  
[www.novusbio.com/761-0015](http://www.novusbio.com/761-0015)

Updated 3/16/2026 v.20.1

**Earn rewards for product  
reviews and publications.**

Submit a publication at [www.novusbio.com/publications](http://www.novusbio.com/publications)

Submit a review at [www.novusbio.com/reviews/destination/761-0015](http://www.novusbio.com/reviews/destination/761-0015)



**761-0015****Lightning-Link (R) PE-Cy5.5 Antibody Labeling Kit**

<b>Product Information</b>	
<b>Unit Size</b>	600 ug
<b>Concentration</b>	Concentration is not relevant for this product. Please see the protocols for proper use of this product.
<b>Storage</b>	Store at -20C.
<b>Conjugate</b>	PE/Cy5.5
<b>Product Description</b>	
<b>Description</b>	<p>Lightning-Link antibody labeling kits enable the direct labeling of antibodies, proteins, peptides or other biomolecules for use in R&amp;D applications, drug discovery and the development of diagnostic kits (See protocol for further information).</p> <p>Our PE/Cy5.5 antibody labeling kit enables the direct conjugation of the PE/Cy5.5 tandem dye to any biomolecule with an available amine group. The researcher simply pipettes the antibody or other biomolecule into the vial of Lightning-Link label and incubates for 3 hours.</p> <p><b>Features</b>Quick and easy to use <b>Benefits</b>Save time, no special knowledge required No separation steps 100% recovery - no antibody/protein loss Can be used in a wide range of applications <b>Flexible</b>Freeze dried Ships at ambient temperature, long shelf-life <b>Fully scalable</b> (10 ug to 1 g or more) <b>Easy transfer</b> from R&amp;D to manufacturing <b>Stringently QC tested</b> <b>Consistent high quality</b>, excellent batch-to-batch reproducibility <b>Large number of labels available</b> <b>Experimental flexibility</b> <b>Reliable</b>: nearly 300 references <b>Successfully used in many fields of research</b></p> <p>PE/Cy5.5 is a tandem conjugate. The PE has three maximal absorbance values of 498, 544 and 566nm with the optimal depending on the application. The PE functions as an energy donor for the Cy5.5. Energy is transferred from the PE to the Cy5.5 via energy resonance transfer. The Cy5.5 emits the energy received from the PE in the form of long wavelength light at 700nm.</p> <p>Learn more about Lightning-Link™ Conjugation Kits by reading <a href="#">FAQs</a></p> <p>For more information please check out these useful links! <a href="#">Antibody Labeling Guide</a> <a href="#">Antibody Conjugation Illustrated Assay</a></p>
<b>Kit Components</b>	1 or 3 glass vial(s) of Lightning-Link mix, 1 vial of LL-Modifier reagent, 1 vial of LL-Quencher reagent
<b>Notes</b>	<p>This product is manufactured by Abcam and distributed by Novus Biologicals.</p> <p>This product is for research use only and is not approved for use in humans or in clinical diagnosis. This product is guaranteed for 1 year from date of receipt and this statement overrides any mentioned guarantee period on the limitations section of this products datasheet. Please contact <a href="mailto:technical@novusbio.com">technical@novusbio.com</a> with questions.</p>
<b>Product Application Details</b>	
<b>Applications</b>	Flow Cytometry
<b>Recommended Dilutions</b>	Flow Cytometry

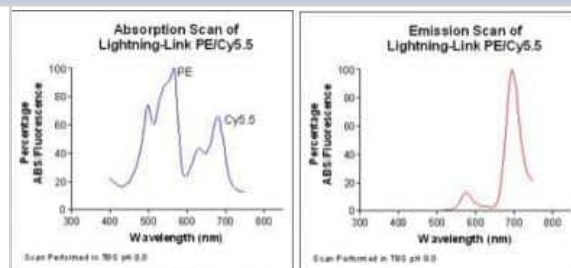


**Application Notes**

By circumventing the desalting or dialysis steps that commonly interrupt traditional antibody conjugation procedures, LightningLink technology can be used to label both small (e.g. 10 ug) and large quantities of primary antibodies with ease. Batch-to-batch variation upon scale up is minimal as the process is so simple, and recoveries are always 100%. This kit can be used to label up to 600 ug of antibody, and is supplied in one vial.

**Images**

Flow Cytometry: Lightning-Link PE-Cy5.5 Antibody Labeling Kit [761-0015]

**Publications**

Wang S, Wang B, Pan N et al. Differentiation of human induced pluripotent stem cells to mature functional Purkinje neurons Sci Rep. 2015-03-18 [PMID: 25782665] (FLOW)

Hesdorffer CS, Malchinkhuu E, Biragyn A et al. Distinctive immunoregulatory effects of adenosine on T cells of older humans The FASEB Journal 2010-01-01 [PMID: 22121051]

Olling A, Goy S, Hoffmann F et al. The Repetitive Oligopeptide Sequences Modulate Cytopathic Potency but Are Not Crucial for Cellular Uptake of Clostridium difficile Toxin A. PLoS One 2011-01-01 [PMID: 21445253] (FLOW)

Carter CC, Onafuwa-Nuga A, McNamara LA et al. HIV-1 Infects Multipotent Progenitor Cells Causing Cell Death and Establishing Latent Cellular Reservoirs Nature Medicine 2010-01-01 [PMID: 20208541]

Al-Dujaili EA, Mullins LJ, Bailey MA, Kenyon CJ. Development of a Highly Sensitive ELISA for Aldosterone in Mouse Urine: Validation in Physiological Pathophysiological States of Aldosterone Excess Depletion - . Steroids. 2008-01-01 [PMID: 19162057]

Bao S, Wu Q, Li Z et al. Targeting Cancer Stem Cells through L1CAM Suppresses Glioma Growth. Cancer Res 2008-01-01 [PMID: 18676824]

Marr AK, Jenssen H, Moniri MR et al. Bovine lactoferrin and lactoferricin interfere with intracellular trafficking of Herpes simplex virus-1. Biochemie 2009-01-01 [PMID: 18573311]

Velappan N, Clements J, Kiss C et al. Fluorescence linked immunosorbant assays using microtiter plates. J Immunol Methods. 2013-01-01 [PMID: 18514691]



### **Novus Biologicals USA**

10730 E. Briarwood Avenue  
Centennial, CO 80112  
USA  
Phone: 303.730.1950  
Toll Free: 1.888.506.6887  
Fax: 303.730.1966  
nb-customerservice@bio-techne.com

### **Bio-Techne Canada**

21 Canmotor Ave  
Toronto, ON M8Z 4E6  
Canada  
Phone: 905.827.6400  
Toll Free: 855.668.8722  
Fax: 905.827.6402  
canada.inquires@bio-techne.com

### **Bio-Techne Ltd**

19 Barton Lane  
Abingdon Science Park  
Abingdon, OX14 3NB, United Kingdom  
Phone: (44) (0) 1235 529449  
Free Phone: 0800 37 34 15  
Fax: (44) (0) 1235 533420  
info.EMEA@bio-techne.com

### **General Contact Information**

www.novusbio.com  
Technical Support: nb-technical@bio-techne.com  
Orders: nb-customerservice@bio-techne.com  
General: novus@novusbio.com

### **Limitations**

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Kits are guaranteed for 6 months from date of receipt.

For more information on our 100% guarantee, please visit [www.novusbio.com/guarantee](http://www.novusbio.com/guarantee)

Earn gift cards/discounts by submitting a review: [www.novusbio.com/reviews/submit/761-0015](http://www.novusbio.com/reviews/submit/761-0015)

Earn gift cards/discounts by submitting a publication using this product:  
[www.novusbio.com/publications](http://www.novusbio.com/publications)

