

Product Datasheet

Lightning-Link (R) PE-Cy7 Antibody Labeling Kit 762-0005

Unit Size: 60 ug

Store at -20C.

www.novusbio.com



technical@novusbio.com

Publications: 14

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:
www.novusbio.com/762-0005

Updated 8/7/2024 v.20.1

**Earn rewards for product
reviews and publications.**

Submit a publication at www.novusbio.com/publications

Submit a review at www.novusbio.com/reviews/destination/762-0005



762-0005**Lightning-Link (R) PE-Cy7 Antibody Labeling Kit**


Product Information	
Unit Size	60 ug
Concentration	Concentration is not relevant for this product. Please see the protocols for proper use of this product.
Storage	Store at -20C.
Conjugate	PE/Cy7
Product Description	
Description	<p>Lightning-Link antibody labeling kits enable the direct labeling of antibodies, proteins, peptides or other biomolecules for use in R&D applications, drug discovery and the development of diagnostic kits (See protocol for further information).</p> <p>Our PE/Cy7 antibody labeling kit enables the direct conjugation of the PE/Cy7 tandem dye to any biomolecule with an available amine group. The researcher simply pipettes the antibody or other biomolecule into the vial of Lightning-Link label and incubates for 3 hours.</p> <p>FeaturesQuick and easy to use BenefitsSave time, no special knowledge required No separation steps 100% recovery - no antibody/protein loss Can be used in a wide range of applications FlexibleFreeze dried Ships at ambient temperature, long shelf-life Fully scalable (10 ug to 1 g or more) Easy transfer from R&D to manufacturing Stringently QC tested Consistent high quality, excellent batch-to-batch reproducibility Large number of labels available Experimental flexibility Reliable: nearly 300 references Successfully used in many fields of research</p> <p>PE/Cy7 is a tandem conjugate. The PE has three maximal absorbance values of 498, 544 and 566nm with the optimal depending on the application. The PE functions as an energy donor for the Cy7. Energy is transferred from the PE to the Cy7 via energy resonance transfer. The Cy7 emits the energy received from the PE in the form of long wavelength light at 782nm.</p> <p>Learn more about Lightning-Link™ Conjugation Kits by reading FAQs</p> <p>For more information please check out these useful links! Antibody Labeling Guide Antibody Conjugation Illustrated Assay</p>
Kit Components	1 or 3 glass vial(s) of Lightning-Link mix, 1 vial of LL-Modifier reagent, 1 vial of LL-Quencher reagent
Notes	<p>This product is manufactured by Abcam and distributed by Novus Biologicals.</p> <p>This product is for research use only and is not approved for use in humans or in clinical diagnosis. This product is guaranteed for 1 year from date of receipt and this statement overrides any mentioned guarantee period on the limitations section of this products datasheet. Please contact technical@novusbio.com with questions.</p>
Product Application Details	
Applications	Flow Cytometry
Recommended Dilutions	Flow Cytometry

Application Notes

By circumventing the desalting or dialysis steps that commonly interrupt traditional antibody conjugation procedures, LightningLink technology can be used to label both small (e.g. 10 ug) and large quantities of primary antibodies with ease. Batch-to-batch variation upon scale up is minimal as the process is so simple, and recoveries are always 100%. This kit can be used to label up to 60 ug of antibody, and is supplied in one vial.

Images

Flow Cytometry: Lightning-Link PE-Cy7 Antibody Labeling Kit [762-0005]

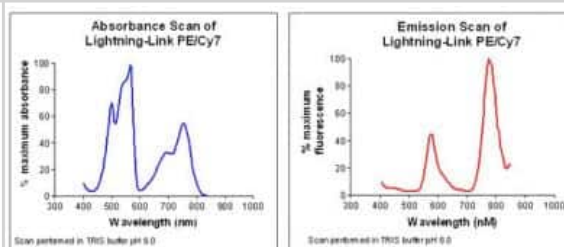

Publications

Yazicioglu YF, Marin E, Sandhu C et al. Dynamic mitochondrial transcription and translation in B cells control germinal center entry and lymphomagenesis *Nature Immunology* 2023-06-01 [PMID: 37095377]

Titus HE, Xu H, Robinson AP et al. Repurposing the cardiac glycoside digoxin to stimulate myelin regeneration in chemically-induced and immune-mediated mouse models of multiple sclerosis *Glia* 2022-07-09 [PMID: 35809238]

Shanmugasundaram R, Acevedo K, Mortada M Et al. Effects of *Salmonella enterica* ser. Enteritidis and Heidelberg on host CD4+CD25+ regulatory T cell suppressive immune responses in chickens *PLOS ONE* 2021-11-29 [PMID: 34843525]

Okagawa T, Konnai S, Nishimori A et al. Bovine immunoinhibitory receptors contribute to the suppression of *Mycobacterium avium* subsp. *paratuberculosis*-specific T-cell responses *Infect Immun.* 2015-10-19 [PMID: 26483406] (FLOW)

Guerra-Perez N, Aravantinou M, Veglia F et al. Rectal HSV-2 Infection May Increase Rectal SIV Acquisition Even in the Context of SIVdeltanef Vaccination *PLoS One* 2016-02-17 [PMID: 26886938] (FLOW)

Boonpiyathad T, Meyer N, Moniuszko M et al. High-dose bee venom exposure induces similar tolerogenic B cell responses in patients and healthy beekeepers *Allergy* 2016-06-24 [PMID: 27341567] (FLOW)

Zhou M, Wang T, Lai H et al. Targeting of the deubiquitinase USP9X attenuates B-cell acute lymphoblastic leukemia cell survival and overcomes glucocorticoid resistance. *Biochem Biophys Res Commun.* 2015-01-01 [PMID: 25735983]

Rodgers JM, Robinson AP, Rosler ES et al. IL-17A activates ERK1/2 and enhances differentiation of oligodendrocyte progenitor cells. *Glia* 2014-01-01 [PMID: 25557204]

Goode D, Truong R, Villegas G et al. Driven Increase in the Expression of $\alpha 4\beta 7$ Correlates with Increased Susceptibility to Vaginal SHIVSF162P3 Infection *PLoS Pathog* 2014-12-18 [PMID: 25521298] (FLOW)

Robinson AP, Rodgers JM, Goings GE, Miller SD. Characterization of Oligodendroglial Populations in Mouse Demyelinating Disease Using Flow Cytometry: Clues for MS Pathogenesis. *PLoS One* 2014-01-01 [PMID: 25247590] (FLOW)

Gallo S, Gatti S, Sala V et al. Agonist antibodies activating the Met receptor protect cardiomyoblasts from cobalt chloride-induced apoptosis and autophagy. *Cell Death Dis* 2014-01-01 [PMID: 24743740] (FLOW)

Hartley AN, Cooley G, Gwyn S et al. Frequency of IFN gamma-producing T cells correlates with seroreactivity and activated T cells during canine *Trypanosoma cruzi* infection. *Vet Res* 2014-01-01 [PMID: 24456537] (FLOW)

More publications at <http://www.novusbio.com/762-0005>





Novus Biologicals USA

10730 E. Briarwood Avenue
Centennial, CO 80112

USA

Phone: 303.730.1950

Toll Free: 1.888.506.6887

Fax: 303.730.1966

nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave

Toronto, ON M8Z 4E6

Canada

Phone: 905.827.6400

Toll Free: 855.668.8722

Fax: 905.827.6402

canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane

Abingdon Science Park

Abingdon, OX14 3NB, United Kingdom

Phone: (44) (0) 1235 529449

Free Phone: 0800 37 34 15

Fax: (44) (0) 1235 533420

info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com

Technical Support: nb-technical@bio-techne.com

Orders: nb-customerservice@bio-techne.com

General: novus@novusbio.com

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Kits are guaranteed for 6 months from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/762-0005

Earn gift cards/discounts by submitting a publication using this product:
www.novusbio.com/publications

