

Product Datasheet

Lightning-Link (R) PE-Cy7 Antibody Labeling Kit 762-0030

Unit Size: 3 x 10ug Reaction

Store at -20C.

www.novusbio.com



technical@novusbio.com

Publications: 10


Protocols, Publications, Related Products, Reviews, Research Tools and Images at:
www.novusbio.com/762-0030

Updated 10/23/2024 v.20.1

**Earn rewards for product
reviews and publications.**

Submit a publication at www.novusbio.com/publications

Submit a review at www.novusbio.com/reviews/destination/762-0030

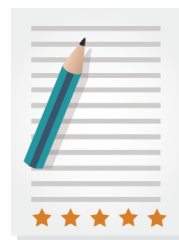

762-0030**Lightning-Link (R) PE-Cy7 Antibody Labeling Kit**

Product Information	
Unit Size	3 x 10ug Reaction
Concentration	Concentration is not relevant for this product. Please see the protocols for proper use of this product.
Storage	Store at -20C.
Conjugate	PE/Cy7
Product Description	
Description	Lightning-Link antibody labeling kits enable the direct labeling of antibodies, proteins, peptides or other biomolecules for use in R&D applications, drug discovery and the development of diagnostic kits (See protocol for further information). Our PE/Cy7 antibody labeling kit enables the direct conjugation of the PE/Cy7 tandem dye to any biomolecule with an available amine group. The researcher simply pipettes the antibody or other biomolecule into the vial of Lightning-Link label and incubates for 3 hours. FeaturesQuick and easy to use BenefitsSave time, no special knowledge required No separation steps 100% recovery - no antibody/protein loss Can be used in a wide range of applications Flexible Freeze dried Ships at ambient temperature, long shelf-life Fully scalable (10 ug to 1 g or more) Easy transfer from R&D to manufacturing Stringently QC tested Consistent high quality, excellent batch-to-batch reproducibility Large number of labels available Experimental flexibility Reliable: nearly 300 references Successfully used in many fields of research PE/Cy7 is a tandem conjugate. The PE has three maximal absorbance values of 498, 544 and 566nm with the optimal depending on the application. The PE functions as an energy donor for the Cy7. Energy is transferred from the PE to the Cy7 via energy resonance transfer. The Cy7 emits the energy received from the PE in the form of long wavelength light at 782nm. Learn more about Lightning-Link™ Conjugation Kits by reading FAQs For more information please check out these useful links! Antibody Labeling Guide Antibody Conjugation Illustrated Assay
Kit Components	1 or 3 glass vial(s) of Lightning-Link mix, 1 vial of LL-Modifier reagent, 1 vial of LL-Quencher reagent
Notes	This product is manufactured by Abcam and distributed by Novus Biologicals. This product is for research use only and is not approved for use in humans or in clinical diagnosis. This product is guaranteed for 1 year from date of receipt and this statement overrides any mentioned guarantee period on the limitations section of this products datasheet. Please contact technical@novusbio.com with questions.
Product Application Details	
Applications	Flow Cytometry
Recommended Dilutions	Flow Cytometry




Application Notes

By circumventing the desalting or dialysis steps that commonly interrupt traditional antibody conjugation procedures, LightningLink technology can be used to label both small (e.g. 10 ug) and large quantities of primary antibodies with ease. Batch-to-batch variation upon scale up is minimal as the process is so simple, and recoveries are always 100%.

Images

Flow Cytometry: Lightning-Link PE-Cy7 Antibody Labeling Kit [762-0030]

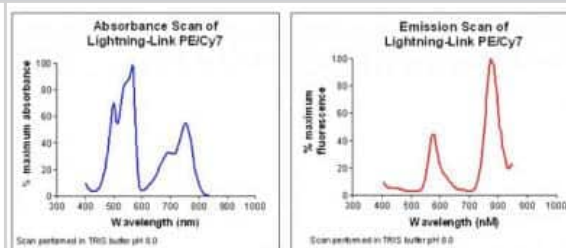**Publications**

Okagawa T, Konnai S, Nishimori A et al. Bovine immunoinhibitory receptors contribute to the suppression of Mycobacterium avium subsp.paratuberculosis-specific T-cell responses Infect Immun. 2015-10-19 [PMID: 26483406] (FLOW)

Guerra-Perez N, Aravantinou M, Veglia F et al. Rectal HSV-2 Infection May Increase Rectal SIV Acquisition Even in the Context of SIVdeltanef Vaccination PLoS One 2016-02-17 [PMID: 26886938] (FLOW)

Boonpiyathad T, Meyer N, Moniuszko M et al. High-dose bee venom exposure induces similar tolerogenic B cell responses in patients and healthy beekeepers Allergy 2016-06-24 [PMID: 27341567] (FLOW)

Zhou M, Wang T, Lai H et al. Targeting of the deubiquitinase USP9X attenuates B-cell acute lymphoblastic leukemia cell survival and overcomes glucocorticoid resistance. Biochem Biophys Res Commun. 2015-01-01 [PMID: 25735983]

Rodgers JM, Robinson AP, Rosler ES et al. IL-17A activates ERK1/2 and enhances differentiation of oligodendrocyte progenitor cells. Glia 2014-01-01 [PMID: 25557204]

Goode D, Truong R, Villegas G et al. Driven Increase in the Expression of a4Beta7 Correlates with Increased Susceptibility to Vaginal SHIVSF162P3 Infection PLoS Pathog 2014-12-18 [PMID: 25521298] (FLOW)

Robinson AP, Rodgers JM, Goings GE, Miller SD. Characterization of Oligodendroglial Populations in Mouse Demyelinating Disease Using Flow Cytometry: Clues for MS Pathogenesis. PLoS One 2014-01-01 [PMID: 25247590] (FLOW)

Gallo S, Gatti S, Sala V et al. Agonist antibodies activating the Met receptor protect cardiomyoblasts from cobalt chloride-induced apoptosis and autophagy. Cell Death Dis 2014-01-01 [PMID: 24743740] (FLOW)

Hartley AN, Cooley G, Gwyn S et al. Frequency of IFN gamma-producing T cells correlates with seroreactivity and activated T cells during canine Trypanosoma cruzi infection. Vet Res 2014-01-01 [PMID: 24456537] (FLOW)

Ikebuchi R, Konnai S, Okagawa T et al. Influence of PD-L1 cross-linking on cell death in PD-L1-expressing cell lines and bovine lymphocytes. Immunology 2013-01-01 [PMID: 24405267]



Novus Biologicals USA

10730 E. Briarwood Avenue
Centennial, CO 80112
USA
Phone: 303.730.1950
Toll Free: 1.888.506.6887
Fax: 303.730.1966
nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave
Toronto, ON M8Z 4E6
Canada
Phone: 905.827.6400
Toll Free: 855.668.8722
Fax: 905.827.6402
canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane
Abingdon Science Park
Abingdon, OX14 3NB, United Kingdom
Phone: (44) (0) 1235 529449
Free Phone: 0800 37 34 15
Fax: (44) (0) 1235 533420
info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com
Technical Support: nb-technical@bio-techne.com
Orders: nb-customerservice@bio-techne.com
General: novus@novusbio.com

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Kits are guaranteed for 6 months from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/762-0030

Earn gift cards/discounts by submitting a publication using this product:
www.novusbio.com/publications

