

Product Datasheet

Lightning-Link (R) PerCP-Cy5.5 Antibody Labeling Kit 763-0030

Unit Size: 3 x 10ug Reaction

Store at -20C.

www.novusbio.com



technical@novusbio.com

Publications: 14

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:
www.novusbio.com/763-0030

Updated 10/23/2024 v.20.1

Earn rewards for product
reviews and publications.

Submit a publication at www.novusbio.com/publications

Submit a review at www.novusbio.com/reviews/destination/763-0030



763-0030**Lightning-Link (R) PerCP-Cy5.5 Antibody Labeling Kit**

Product Information	
Unit Size	3 x 10ug Reaction
Concentration	Concentration is not relevant for this product. Please see the protocols for proper use of this product.
Storage	Store at -20C.
Conjugate	PerCP/Cy5.5
Product Description	
Description	Lightning-Link antibody labeling kits enable the direct labeling of antibodies, proteins, peptides or other biomolecules for use in R&D applications, drug discovery and the development of diagnostic kits (See protocol for further information). Our PerCP/Cy5.5 antibody labeling kit enables the direct conjugation of the PerCP/Cy5.5 tandem dye to any biomolecule with an available amine group. The researcher simply pipettes the antibody or other biomolecule into the vial of Lightning-Link label and incubates for 3 hours. FeaturesQuick and easy to use BenefitsSave time, no special knowledge required No separation steps 100% recovery - no antibody/protein loss Can be used in a wide range of applications Flexible Freeze dried Ships at ambient temperature, long shelf-life Fully scalable (10 ug to 1 g or more) Easy transfer from R&D to manufacturing Stringently QC tested Consistent high quality, excellent batch-to-batch reproducibility Large number of labels available Experimental flexibility Reliable: nearly 300 references Successfully used in many fields of research PerCP/Cy5.5 is a tandem label. The PerCP excited at 484nm and functions as an energy donor for the Cy5.5. Energy is transferred from the PerCP to the Cy5.5 via energy resonance transfer. The Cy5.5 emits the energy received from the PerCP in the form of long wavelength light at 692nm. Learn more about Lightning-Link™ Conjugation Kits by reading FAQs For more information please check out these useful links! Antibody Labeling Guide Antibody Conjugation Illustrated Assay
Kit Components	1 or 3 glass vial(s) of Lightning-Link mix, 1 vial of LL-Modifier reagent, 1 vial of LL-Quencher reagent
Notes	This product is manufactured by Abcam and distributed by Novus Biologicals. This product is for research use only and is not approved for use in humans or in clinical diagnosis. This product is guaranteed for 1 year from date of receipt and this statement overrides any mentioned guarantee period on the limitations section of this products datasheet. Please contact technical@novusbio.com with questions.
Product Application Details	
Applications	Electron Microscopy, Flow Cytometry, Flow (Cell Surface)
Recommended Dilutions	Flow Cytometry, Electron Microscopy, Flow (Cell Surface)




Application Notes

PerCP/Cy5.5 is a tandem conjugate. In our Lightning-Link PerCP/Cy5.5 the best result is obtained when the tandem is excited at 488nm, wavelength at which the PerCP has its maximum absorption. The excited PerCP functions as an energy donor, where the energy is transferred via FRET (Fluorescence Resonance Energy Transfer) to the Cy5.5. The Cy5.5 emits then the energy received from the PE in the form of long wavelength light at 700nm. Use in FLOW cytometry reported in scientific literature (PMID 26483406). This kit is supplied with 3 vials, each suitable for labeling up to 10 ug of antibody.

Images

Flow Cytometry: Lightning-Link PerCP-Cy5.5 Antibody Labeling Kit [763-0030] - Mouse anti-human CD3 was conjugated with PerCP/Cy5.5 using Lightning-Link kit. The conjugated antibody was then used to stain human peripheral blood lymphocytes, followed by analysis with flow cytometry. (Blue line - negative control; red line - positive staining).

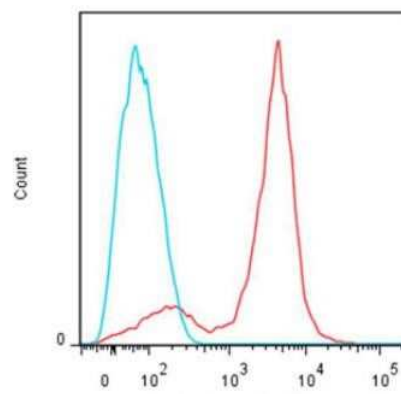

Publications

Goto S, Konnai S, Okagawa T et al. Increase of cells expressing PD-1 and PD-L1 and enhancement of IFN- γ production via PD-1/PD-L1 blockade in bovine mycoplasmosis. *Immun Inflamm Dis*. 2017-01-01 [PMID: 28544524]

Kjaerup RB, Juul-Madsen HR, Norup LR et al. Comparison of growth performance and immune parameters of three commercial chicken lines used in organic production. *Vet Immunol Immunopathol*. 2017-01-01 [PMID: 28494932]

Nishimori A, Konnai S, Okagawa T et al. In vitro and in vivo antiviral activity of an anti-programmed death-ligand 1 (PD-L1) rat-bovine chimeric antibody against bovine leukemia virus infection. *PLoS One*. 2017-04-26 [PMID: 28445479]

Okagawa T, Konnai S, Deringer JR et al. Cooperation of PD-1 and LAG-3 contributes to T-cell exhaustion in *Anaplasma marginale*-infected cattle *Infect Immun*. 2016-07-18 [PMID: 27430272] (FLOW)

Kitazawa Y, Ueta H, Hunig T et al. A novel multicolor immunostaining method using ethynyl deoxyuridine for analysis of in situ immunoproliferative response *Histochem Cell Biol*. 2015-09-01 [PMID: 25976155] (FLOW)

Vremec D. The Isolation and Enrichment of Large Numbers of Highly Purified Mouse Spleen Dendritic Cell Populations and Their In Vitro Equivalents *Methods Mol Biol*. 2016-01-01 [PMID: 27142009] (FLOW)

Riley SP, Fish AI, Garza DA et al. Non-selective Persistence of a *Rickettsia conorii* Extrachromosomal Plasmid During Mammalian Infection *Infect Immun*. 2016-01-11 [PMID: 26755154] (FLOW)

Dooley J, Garcia-Perez JE, Sreenivasan J et al. The microRNA-29 Family Dictates the Balance Between Homeostatic and Pathological Glucose Handling in Diabetes and Obesity *Diabetes* 2016-01-01 [PMID: 26696639] (FLOW)

Okagawa T, Konnai S, Nishimori A et al. Bovine immunoinhibitory receptors contribute to the suppression of *Mycobacterium avium* subsp. *paratuberculosis*-specific T-cell responses *Infect Immun*. 2015-10-19 [PMID: 26483406] (FLOW)

Goode D, Truong R, Villegas G et al. Driven Increase in the Expression of $\alpha 4\beta 7$ Correlates with Increased Susceptibility to Vaginal SHIVSF162P3 Infection *PLoS Pathog* 2014-12-18 [PMID: 25521298] (FLOW)

Slota C, Shi A, Chen G et al. Norepinephrine preferentially modulates memory CD8 T cell function inducing inflammatory cytokine production and reducing proliferation in response to activation. *Brain Behav Immun*. 2015-01-01 [PMID: 25653192] (FLOW)

Robinson AP, Rodgers JM, Goings GE, Miller SD. Characterization of Oligodendroglial Populations in Mouse Demyelinating Disease Using Flow Cytometry: Clues for MS Pathogenesis. *PLoS One* 2014-01-01 [PMID: 25247590] (FLOW)

More publications at <http://www.novusbio.com/763-0030>





Novus Biologicals USA

10730 E. Briarwood Avenue
Centennial, CO 80112
USA
Phone: 303.730.1950
Toll Free: 1.888.506.6887
Fax: 303.730.1966
nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave
Toronto, ON M8Z 4E6
Canada
Phone: 905.827.6400
Toll Free: 855.668.8722
Fax: 905.827.6402
canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane
Abingdon Science Park
Abingdon, OX14 3NB, United Kingdom
Phone: (44) (0) 1235 529449
Free Phone: 0800 37 34 15
Fax: (44) (0) 1235 533420
info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com
Technical Support: nb-technical@bio-techne.com
Orders: nb-customerservice@bio-techne.com
General: novus@novusbio.com

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Kits are guaranteed for 6 months from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/763-0030

Earn gift cards/discounts by submitting a publication using this product:
www.novusbio.com/publications

