

Product Datasheet

Collagen I alpha 1 Antibody - BSA Free NBP1-77457

Unit Size: 0.1 ml

Aliquot and store at -20C or -80C. Avoid freeze-thaw cycles.

www.novusbio.com



technical@novusbio.com

Publications: 3

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:
www.novusbio.com/NBP1-77457

Updated 9/9/2025 v.20.1

Earn rewards for product
reviews and publications.

Submit a publication at www.novusbio.com/publications

Submit a review at www.novusbio.com/reviews/destination/NBP1-77457



NBP1-77457

Collagen I alpha 1 Antibody - BSA Free

Product Information	
Unit Size	0.1 ml
Concentration	1.0 mg/ml
Storage	Aliquot and store at -20C or -80C. Avoid freeze-thaw cycles.
Clonality	Polyclonal
Preservative	0.05% Sodium Azide
Isotype	IgG
Purity	Immunogen affinity purified
Buffer	PBS and 30% Glycerol
Target Molecular Weight	139 kDa


Product Description	
Description	Novus Biologicals Rabbit Collagen I alpha 1 Antibody - BSA Free (NBP1-77457) is a polyclonal antibody validated for use in IHC and ICC/IF. Anti-Collagen I alpha 1 Antibody: Cited in 2 publications. All Novus Biologicals antibodies are covered by our 100% guarantee.
Host	Rabbit
Gene ID	1277
Gene Symbol	COL1A1
Species	Human, Mouse
Immunogen	This Collagen I alpha 1 antibody was raised against a synthetic peptide made to an internal portion of the human Collagen I protein (between residues 150-200) [UniProt# P02452]

Product Application Details	
Applications	Immunohistochemistry-Paraffin, Immunocytochemistry/ Immunofluorescence, Immunohistochemistry, Western Blot (Negative)
Recommended Dilutions	Immunohistochemistry 1:200, Immunocytochemistry/ Immunofluorescence 10-20 ug/ml, Immunohistochemistry-Paraffin 1:200, Western Blot (Negative)
Application Notes	This Collagen 1 antibody is useful for IHC-paraffin embedded sections.




Images

Immunohistochemistry: Collagen I alpha 1 Antibody [NBP1-77457] - Collagen I Antibody [NBP1-77457] - Analysis of Collagen 1 in mouse epidermis using DAB with hematoxylin counterstain.

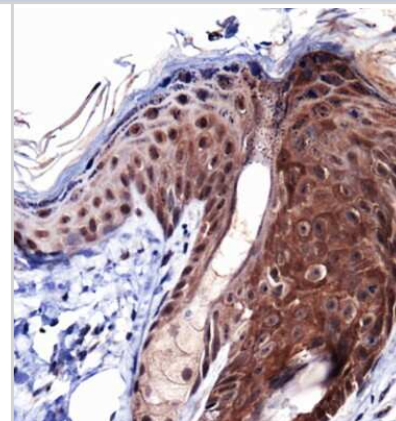

Immunocytochemistry/Immunofluorescence: Collagen I alpha 1 Antibody [NBP1-77457] - Collagen 1 antibody was tested in Hela cells with DyLight 488 (green). Nuclei and alpha-tubulin were counterstained with DAPI (blue) and DyLight 550 (red).

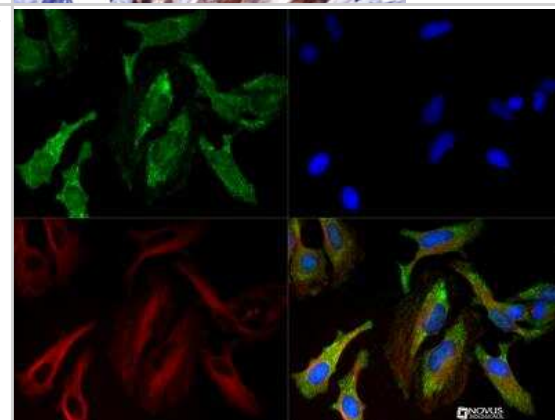

Publications

Lee Y, Lee J, Han H et al. Gingival biotype modification with collagen matrix or autogenous subepithelial connective tissue graft: Histologic and volumetric analyses in a beagle model Heliyon 2023-04-01 [PMID: 37151660]

Li M, Zhou S, Wu Y et al. Prenatal exposure to propylparaben at human-relevant doses accelerates ovarian aging in adult mice Environmental Pollution 2021-09-01 [PMID: 33957517]

S. Kara, M. Akbaba, S. Kul, K. Bakır Is it possible to make early wound age estimation by immunohistochemical methods? Rom J Leg Med 2016-01-01 (IF/IHC, Mouse)

Procedures

Immunohistochemistry-Paraffin Embedded Sections protocol specific for Collagen I antibody (NBP1-77457)

Collagen I alpha 1 Antibody:

Immunohistochemistry-Paraffin Embedded Sections

Antigen Unmasking:

Bring slides to a boil in 10 mM sodium citrate buffer (pH 6.0) then maintain at a sub-boiling temperature for 10 minutes. Cool slides on bench-top for 30 minutes.

Staining:

1. Wash sections in deionized water three times for 5 minutes each.
2. Wash sections in wash buffer for 5 minutes.
3. Block each section with 100-400 ul blocking solution for 1 hour at room temperature.
4. Remove blocking solution and add 100-400 ul diluted primary antibody. Incubate overnight at 4C.
5. Remove antibody solution and wash sections in wash buffer three times for 5 minutes each.
6. Add 100-400 ul biotinylated diluted secondary antibody. Incubate 30 minutes at room temperature.
7. Remove secondary antibody solution and wash sections three times with wash buffer for 5 minutes each.
8. Add 100-400 ul Streptavidin-HRP reagent to each section and incubate for 30 minutes at room temperature.
9. Wash sections three times in wash buffer for 5 minutes each.
10. Add 100-400 ul DAB substrate to each section and monitor staining closely.
11. As soon as the sections develop, immerse slides in deionized water.
12. Counterstain sections in hematoxylin.
13. Wash sections in deionized water two times for 5 minutes each.
14. Dehydrate sections.
15. Mount coverslips.

*The above information is only intended as a guide. The researcher should determine what protocol best meets their needs. Please follow safe laboratory procedures.





Novus Biologicals USA

10730 E. Briarwood Avenue
Centennial, CO 80112
USA
Phone: 303.730.1950
Toll Free: 1.888.506.6887
Fax: 303.730.1966
nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave
Toronto, ON M8Z 4E6
Canada
Phone: 905.827.6400
Toll Free: 855.668.8722
Fax: 905.827.6402
canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane
Abingdon Science Park
Abingdon, OX14 3NB, United Kingdom
Phone: (44) (0) 1235 529449
Free Phone: 0800 37 34 15
Fax: (44) (0) 1235 533420
info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com
Technical Support: nb-technical@bio-techne.com
Orders: nb-customerservice@bio-techne.com
General: novus@novusbio.com

Products Related to NBP1-77457

NBP2-33376H	Blue Marker Antibody (6F4-F6) [HRP]
HAF008	Goat anti-Rabbit IgG Secondary Antibody [HRP]
NB7160	Goat anti-Rabbit IgG (H+L) Secondary Antibody [HRP]
NBP2-24891	Rabbit IgG Isotype Control

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NBP1-77457

Earn gift cards/discounts by submitting a publication using this product:
www.novusbio.com/publications

