

# Product Datasheet

## BST2 Antibody (4F6) - BSA Free NBP2-29622

Unit Size: 0.1 mg

Store at 4C short term. Aliquot and store at -20C long term. Avoid freeze-thaw cycles.

[www.novusbio.com](http://www.novusbio.com)



[technical@novusbio.com](mailto:technical@novusbio.com)

### Publications: 1

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:  
[www.novusbio.com/NBP2-29622](http://www.novusbio.com/NBP2-29622)

Updated 9/9/2025 v.20.1

Earn rewards for product  
reviews and publications.

Submit a publication at [www.novusbio.com/publications](http://www.novusbio.com/publications)

Submit a review at [www.novusbio.com/reviews/destination/NBP2-29622](http://www.novusbio.com/reviews/destination/NBP2-29622)



**NBP2-29622**


BST2 Antibody (4F6) - BSA Free

<b>Product Information</b>	
<b>Unit Size</b>	0.1 mg
<b>Concentration</b>	1.0 mg/ml
<b>Storage</b>	Store at 4C short term. Aliquot and store at -20C long term. Avoid freeze-thaw cycles.
<b>Clonality</b>	Monoclonal
<b>Clone</b>	4F6
<b>Preservative</b>	0.05% Sodium Azide
<b>Isotype</b>	IgG1 Kappa
<b>Purity</b>	Protein A or G purified
<b>Buffer</b>	PBS
<b>Product Description</b>	
<b>Description</b>	Novus Biologicals Mouse BST2 Antibody (4F6) - BSA Free (NBP2-29622) is a monoclonal antibody validated for use in IHC, WB and Flow. Anti-BST2 Antibody: Cited in 1 publication. All Novus Biologicals antibodies are covered by our 100% guarantee.
<b>Host</b>	Mouse
<b>Gene ID</b>	684
<b>Gene Symbol</b>	BST2
<b>Species</b>	Human
<b>Immunogen</b>	Full length recombinant human BST2 protein [NCBI NP_004326].
<b>Product Application Details</b>	
<b>Applications</b>	Western Blot, Immunohistochemistry-Paraffin, Flow Cytometry, Flow (Cell Surface), Immunohistochemistry, CyTOF-ready
<b>Recommended Dilutions</b>	Western Blot 1.0 ug/mL, Flow Cytometry 0.5 ug/mL, Immunohistochemistry 5 ug/mL, Immunohistochemistry-Paraffin 5 ug/mL, Flow (Cell Surface), CyTOF-ready
<b>Application Notes</b>	The unprocessed form of BST2 protein is 180 amino acids long (161 amino acids in the cleaved BST) and its cleaved form undergoes further modifications namely disulfide bond formation, GPI-anchor/GPI-like-anchor binding, glycosylation, isopeptide bond formation, lipidation and Ubl conjugation. The predicted molecular weight of unmodified BST2 protein is ~17-18 kDa, whereas the post-translationally modified protein has been documented to run at ~25-36 kDa (glycosylated), 50-55 kDa (ubiquitinated) and 65-70 kDa (dimers). BST2 protein localizes mainly to the cell membranes from where it shuttles to the trans-Golgi network, late endosome and the cytoplasm. In immuno-staining assays, BST2 may develop a membrane-cytoplasmic staining pattern. This antibody is CyTOF ready.




## Images


Western Blot: BST2 Antibody (4F6) [NBP2-29622] - WB analysis of BST2 protein in (A) human heart lysate (B) human ovary lysate and on (C) partial recombinant BST2 protein with monoclonal BST2 antibody (clone 4F6) at a concentration of 1 ug/mL. In the tested lysates, two bands were observed which represents the glycosylated (lower intense bands) and the ubiquitinated (higher weak bands) forms of BST2 protein.

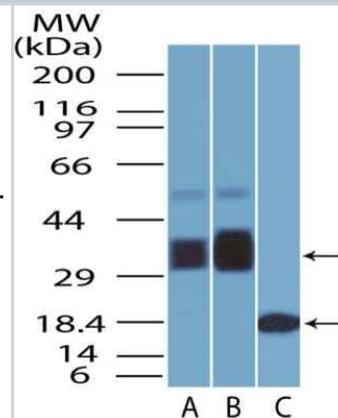


Immunohistochemistry-Paraffin: BST2 Antibody (4F6) [NBP2-29622] - Analysis of BST2 protein in a section of human lung cancer (squamous cell carcinoma) using 5 ug/mL concentration of BST2 antibody (clone 4F6). The carcinoma cells depicted strong membrane-cytoplasmic BST2 positivity with more intense staining in the cellular membranes, and a relatively less intense staining pattern was observed in the tumor stroma.

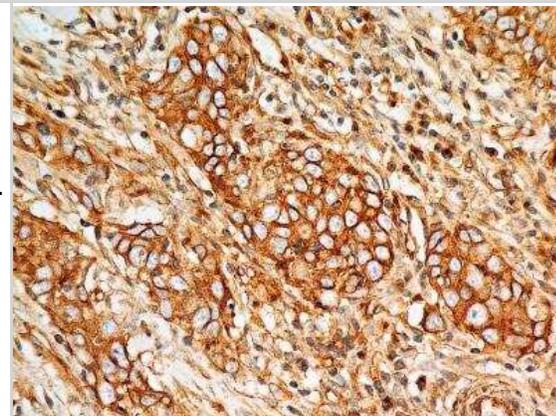

Flow (Cell Surface): BST2 Antibody (4F6) [NBP2-29622] - Flow analysis of Molt4 cells using the BST2 antibody [4F6]. The shaded purple region represents cells alone, green represents the isotype control and pink represents the BST2 antibody.

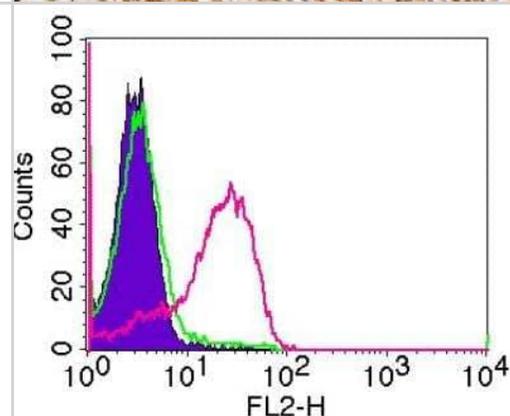

Immunohistochemistry-Paraffin: BST2 Antibody (4F6) [NBP2-29622] - Analysis of BST2 protein in a section of human endometrial carcinoma using 5 ug/mL concentration of BST2 antibody (clone 4F6). The carcinoma cells depicted distinct membrane-cytoplasmic BST2 positivity with more intense staining in the cellular membranes.

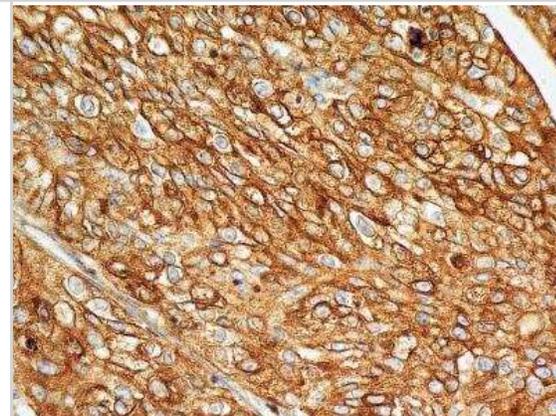

## Publications

Salem DP, Bortolin LT, Gusenleitner D et al. Colocalization of cancer-associated biomarkers on single extracellular vesicles for early detection of cancer. The Journal of molecular diagnostics : JMD 2024-09-24 [PMID: 39326670]

## Procedures

### Western Blot protocol for BST2 Antibody (NBP2-29622)

BST2 Antibody (4F6):

Western Blot Protocol

1. Perform SDS-PAGE on samples to be analyzed, loading 25 ug of total protein per lane.
2. Transfer proteins to membrane according to the instructions provided by the manufacturer of the membrane and transfer apparatus.
3. Stain according to standard Ponceau S procedure (or similar product) to assess transfer success, and mark molecular weight standards where appropriate.
4. Rinse the blot.
5. Block the membrane using standard blocking buffer for at least 1 hour.
6. Wash the membrane in wash buffer three times for 10 minutes each.
7. Dilute anti-BST2 primary antibody in blocking buffer and incubate 1 hour at room temperature.
8. Wash the membrane in wash buffer three times for 10 minutes each.
9. Apply the diluted HRP conjugated secondary antibody in blocking buffer (as per manufacturers instructions) and incubate 1 hour at room temperature.
10. Wash the blot in wash buffer three times for 10 minutes each (this step can be repeated as required to reduce background).
11. Apply the detection reagent of choice in accordance with the manufacturers instructions.

Note: Tween-20 can be added to the blocking or antibody dilution buffer at a final concentration of 0.05-0.2%.



**Immunohistochemistry-Paraffin protocol for BST2 Antibody (NBP2-29622)****BST2 Antibody (4F6):**

1. Deparaffinize the tissue sections by immersing the slides in Xylene with two changes for 10 min each. Sections should not get dried at any stage from this point.
2. Rehydrate the tissue sections by immersing the slides in decreasing grades of ethanol as follows:
  - a. Immerse in 100% ethanol with 2 changes for 5 minutes each
  - b. Immerse in 95% ethanol with 2 changes for 5 minutes each
  - c. Immerse in 90% ethanol for 5 minutes
  - d. Immerse in 70% ethanol for 5 minutes
  - e. Immerse in 50% ethanol for 5 minutes
  - f. Immerse in distilled water for 5 minutes
3. Antigen Retrieval (Microwave Method):
  - a. Immerse the slides in a microwave compatible tray containing 10 mM Sodium Citrate buffer (pH 6.0) with 0.05% Tween 20.
  - b. Boil the slides and maintain the sub-boiling temperature for 5 minutes in the microwave. Thereafter, take out the tray very carefully and cool it at room temperature (RT) for about 30 minutes.
  - c. Wash the slides 3 times, 3 minutes each by immersing them in TBST (Tris Buffered Saline having 0.05% Tween 20).
4. Quenching of Endogenous Peroxidase:
  - a. Incubate the slides in 3% hydrogen peroxide prepared in methanol for 15 minutes (at RT, in dark conditions).
  - b. Wash the slides in TBST 3 times, 3 minutes each.
5. Protein Blocking:
  - a. Incubate the sections with background sniper solution at RT for 15 minutes (Biocare Medicals, USA).
  - b. Wash the sections 3 times, 3 min each by immersing the slides in TBST.
6. Primary Antibody:
  - a. Dilute the primary antibody at 5ug/ml concentration using PBS as a diluent.
  - b. Incubate the sections with diluted primary antibody for 90 minutes at RT in a humidified chamber.
  - c. Thereafter, wash the slides 4 times, 5 minutes each with TBST.
7. Probe (Secondary Reagent):
  - a. Incubate with MACH 1 Mouse probe for 15 minutes at RT.
  - b. Incubate for 30 min at room temperature with HRP-Polymer (Biocare Medical, USA).
  - c. Wash the slides with TBST 4 times, 5 minutes each
8. Chromogen:
  - a. Mix 32ul of DAB Chromogen with 1 ml of DAB substrate buffer (Biocare Medical, USA).
  - b. Apply 200ul DAB mixture/section and incubate at RT in dark conditions (few seconds - 5 minutes).
  - c. As soon as an appropriate color develops, rinse the slides with deionized water (2-3 brief rinses).
9. Counter stain:
  - a. Counter stain with Hematoxylin for 30 seconds (Vector Labs, USA).
  - b. Wash in deionized water for 1-2 minutes to clear the extra stain.
  - c. Incubate the slides in bluing solution or Scott's water twice for 2 minutes each time.
10. Dehydrate the sections in increasing grades of alcohols:
  - a. 50% alcohol for 1 minute
  - b. 70% for 1 minute
  - c. 90% for 1 minute
  - d. 95% for 1 minute
  - e. 100% for 1 minute
  - f. Xylene with 2 changes for 2 minutes each
11. Mount with DPX mount and cover-slip glass (Fisher Scientific, USA), carefully not allowing any air bubbles to enter.

**NOTE:-** This protocol is provided as a reference tool only. Depending upon the type of tissues /tissue processing and reagents employed, the end user will need to optimize the final conditions for achieving an expected staining.



### **Novus Biologicals USA**

10730 E. Briarwood Avenue  
Centennial, CO 80112  
USA  
Phone: 303.730.1950  
Toll Free: 1.888.506.6887  
Fax: 303.730.1966  
nb-customerservice@bio-techne.com

### **Bio-Techne Canada**

21 Canmotor Ave  
Toronto, ON M8Z 4E6  
Canada  
Phone: 905.827.6400  
Toll Free: 855.668.8722  
Fax: 905.827.6402  
canada.inquires@bio-techne.com

### **Bio-Techne Ltd**

19 Barton Lane  
Abingdon Science Park  
Abingdon, OX14 3NB, United Kingdom  
Phone: (44) (0) 1235 529449  
Free Phone: 0800 37 34 15  
Fax: (44) (0) 1235 533420  
info.EMEA@bio-techne.com

### **General Contact Information**

www.novusbio.com  
Technical Support: nb-technical@bio-techne.com  
Orders: nb-customerservice@bio-techne.com  
General: novus@novusbio.com

### **Products Related to NBP2-29622**

---

NBP2-33376H	Blue Marker Antibody (6F4-F6) [HRP]
HAF007	Goat anti-Mouse IgG Secondary Antibody [HRP]
NB7539	Goat anti-Mouse IgG (H+L) Secondary Antibody [HRP]
NBP1-43319-0.5mg	Mouse IgG1 Kappa Isotype Control (P3.6.2.8.1)

---

### **Limitations**

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit [www.novusbio.com/guarantee](http://www.novusbio.com/guarantee)

Earn gift cards/discounts by submitting a review: [www.novusbio.com/reviews/submit/NBP2-29622](http://www.novusbio.com/reviews/submit/NBP2-29622)

Earn gift cards/discounts by submitting a publication using this product:  
[www.novusbio.com/publications](http://www.novusbio.com/publications)

