

# Product Datasheet

## HNF-3 alpha/FoxA1 Antibody (3B11NB) - BSA Free NBP2-45269

Unit Size: 0.1 mg

Store at 4C short term. Aliquot and store at -20C long term. Avoid freeze-thaw cycles.

[www.novusbio.com](http://www.novusbio.com)



[technical@novusbio.com](mailto:technical@novusbio.com)

### Publications: 2

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:  
[www.novusbio.com/NBP2-45269](http://www.novusbio.com/NBP2-45269)

Updated 9/9/2025 v.20.1

Earn rewards for product  
reviews and publications.

Submit a publication at [www.novusbio.com/publications](http://www.novusbio.com/publications)

Submit a review at [www.novusbio.com/reviews/destination/NBP2-45269](http://www.novusbio.com/reviews/destination/NBP2-45269)



**NBP2-45269**

HNF-3 alpha/FoxA1 Antibody (3B11NB) - BSA Free


Product Information	
Unit Size	0.1 mg
Concentration	1.0 mg/ml
Storage	Store at 4C short term. Aliquot and store at -20C long term. Avoid freeze-thaw cycles.
Clonality	Monoclonal
Clone	3B11NB
Preservative	0.02% Sodium Azide
Isotype	IgG2b Kappa
Purity	Protein G purified
Buffer	PBS

Product Description	
Description	Novus Biologicals Mouse HNF-3 alpha/FoxA1 Antibody (3B11NB) - BSA Free (NBP2-45269) is a monoclonal antibody validated for use in IHC, WB, Flow and ICC/IF. Anti-HNF-3 alpha/FoxA1 Antibody: Cited in 2 publications. All Novus Biologicals antibodies are covered by our 100% guarantee.
Host	Mouse
Gene ID	3169
Gene Symbol	FOXA1
Species	Human
Immunogen	Partial recombinant human FOXA1 protein (between amino acids 350-500) [Uniprot: P55317].


Product Application Details	
Applications	Western Blot, Immunohistochemistry-Paraffin, Flow Cytometry, Immunocytochemistry/ Immunofluorescence, Immunohistochemistry
Recommended Dilutions	Western Blot 2-4 ug/ml, Flow Cytometry 25 ug/ml, Immunohistochemistry 5 ug/ml, Immunocytochemistry/ Immunofluorescence 1:25 - 1:50, Immunohistochemistry-Paraffin 5 ug/ml

**Images**

Western Blot: HNF-3 alpha/FoxA1 Antibody (3B11NB) [NBP2-45269] - Whole cell protein from PC3 (1) and MCF7 (2) was separated on a 7.5% gel by SDS-PAGE, transferred to PVDF membrane and blocked in 5% non-fat milk in TBST. The membrane was probed with 4 ug/ml anti-hFOXA1 (3B11NB) in 1% milk and detected with an anti-mouse HRP secondary antibody using chemiluminescence.

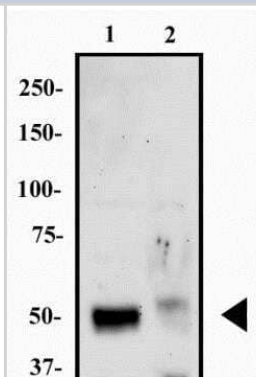


Immunocytochemistry/Immunofluorescence: HNF-3 alpha/FoxA1 Antibody (3B11NB) [NBP2-45269] - FOXA1 (3B11NB) antibody was tested in HepG2 cells at a 1:25 dilution against DyLight 488 (Green) using a 40X objective. Actin and nuclei were counterstained against Phalloidin 568 (Red) and DAPI (Blue), respectively

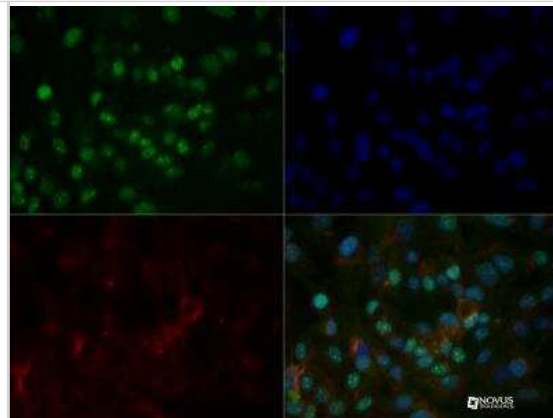


Immunohistochemistry-Paraffin: HNF-3 alpha/FoxA1 Antibody (3B11NB) [NBP2-45269] - Analysis of FOXA1 (3B11NB) in human prostate.




Flow Cytometry: HNF-3 alpha/FoxA1 Antibody (3B11NB) [NBP2-45269] - Analysis using the Azide Free version of NBP2-45269. Staining of LnCAP cells ( $1 \times 10^6$  cells/ml) with FOXA1 antibody, clone 3B11NB (orange), shown with the negative control (blue).

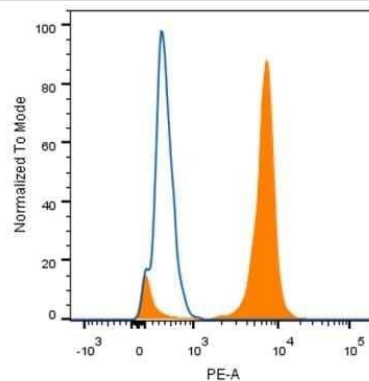


Western Blot: HNF-3 alpha/FoxA1 Antibody (3B11NB) [NBP2-45269] - Analysis of FOXA1 (3B11NB) in partial recombinant FOXA1 protein and PC3 lysate, respectively.

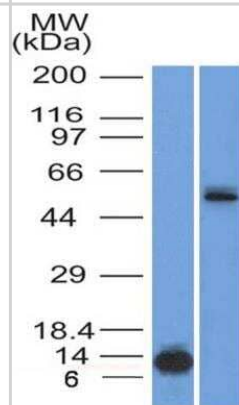


Immunohistochemistry-Paraffin: HNF-3 alpha/FoxA1 Antibody (3B11NB) [NBP2-45269] - Analysis of FOXA1 (3B11NB) in human breast adenocarcinoma.

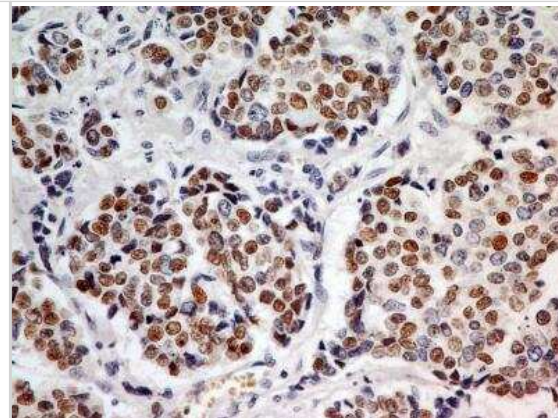

Flow Cytometry: HNF-3 alpha/FoxA1 Antibody (3B11NB) [NBP2-45269] - Analysis using the Azide Free version of NBP2-45269. Staining of HepG2 cells ( $1 \times 10^6$  cells/ml) with FOXA1 antibody, clone 3B11NB (orange), shown with the negative control (blue).

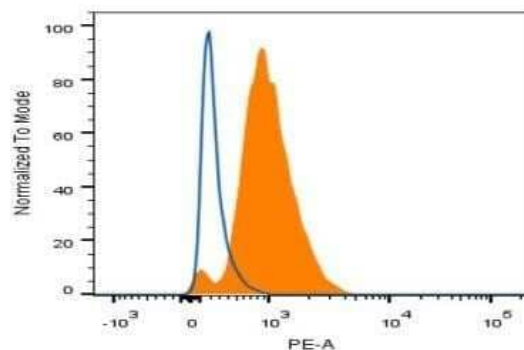

Western Blot: Mouse Monoclonal HNF-3 alpha/FoxA1 Antibody (3B11NB) [NBP2-45269] - 40 ug protein from OUMS29 cell lysates were loaded into gels for western blot experiment. The PVDF membranes were probed with this antibody at 1ug/ml concentration. Bands for HNF-3 alpha/FoxA1 appeared around 53kDa. Image from a verified customer review.

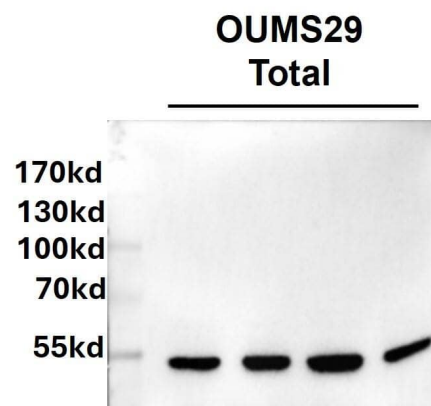

## Publications

Han M, Li F, Zhang Y et al. FOXA2 drives lineage plasticity and KIT pathway activation in neuroendocrine prostate cancer *Cancer cell* 2022-11-14 [PMID: 36332622] (IHC-P, Mouse)

Details:

Dilution used for IHC 1:200

Chen Y L, Hu C M et al. Cellular 5-hydroxymethylcytosine content determines tumorigenic potential and prognosis of pancreatic ductal adenocarcinoma. *Am J Cancer Res* 2019-01-22 [PMID: 30662811] (IF/IHC, Human)

## Procedures

### Western Blot protocol for HNF-3 alpha/FoxA1 Antibody (NBP2-45269)

HNF-3 alpha/FoxA1 Antibody (3B11NB):

1. Perform SDS-PAGE on samples to be analyzed, loading 25 ug of total protein per lane.
2. Transfer proteins to membrane according to the instructions provided by the manufacturer of the membrane and transfer apparatus.
3. Stain according to standard Ponceau S procedure (or similar product) to assess transfer success, and mark molecular weight standards where appropriate.
4. Rinse the blot.
5. Block the membrane using standard blocking buffer for at least 1 hour.
6. Wash the membrane in wash buffer three times for 10 minutes each.
7. Dilute anti-FOXA1 primary antibody in blocking buffer and incubate 1 hour at room temperature.
8. Wash the membrane in wash buffer three times for 10 minutes each.
9. Apply the diluted HRP conjugated secondary antibody in blocking buffer (as per manufacturers instructions) and incubate 1 hour at room temperature.
10. Wash the blot in wash buffer three times for 10 minutes each (this step can be repeated as required to reduce background).
11. Apply the detection reagent of choice in accordance with the manufacturers instructions.

Note: Tween-20 can be added to the blocking or antibody dilution buffer at a final concentration of 0.05-0.2%.

### Immunocytochemistry/Immunofluorescence protocol for HNF-3 alpha/FoxA1 Antibody (NBP2-45269)

HNF-3 alpha/FoxA1 Antibody (3B11NB):

Immunocytochemistry Protocol

Culture cells to appropriate density in 35 mm culture dishes or 6-well plates.

1. Remove culture medium and add 10% formalin to the dish. Fix at room temperature for 30 minutes.
2. Remove the formalin and add ice cold methanol. Incubate for 5-10 minutes.
3. Remove methanol and add washing solution (i.e. PBS). Be sure to not let the specimen dry out. Wash three times for 10 minutes.
4. To block nonspecific antibody binding incubate in 10% normal goat serum from 1 hour to overnight at room temperature.
5. Add primary antibody at appropriate dilution and incubate at room temperature from 2 hours to overnight at room temperature.
6. Remove primary antibody and replace with washing solution. Wash three times for 10 minutes.
7. Add secondary antibody at appropriate dilution. Incubate for 1 hour at room temperature.
8. Remove antibody and replace with wash solution, then wash for 10 minutes. Add Hoechst 33258 to wash solution at 1:25,000 and incubate for 10 minutes. Wash a third time for 10 minutes.
9. Cells can be viewed directly after washing. The plates can also be stored in PBS containing Azide covered in Parafilm (TM). Cells can also be cover-slipped using Fluoromount, with appropriate sealing.

\*The above information is only intended as a guide. The researcher should determine what protocol best meets their needs. Please follow safe laboratory

**Immunohistochemistry-Paraffin protocol for HNF-3 alpha/FoxA1 Antibody (NBP2-45269)****HNF-3 alpha/FoxA1 Antibody (3B11NB):**

1. Deparaffinize the tissue sections by immersing the slides in Xylene with two changes for 10 min each. Sections should not get dried at any stage from this point.
2. Rehydrate the tissue sections by immersing the slides in decreasing grades of ethanol as follows:
  - a. Immerse in 100% ethanol with 2 changes for 5 minutes each
  - b. Immerse in 95% ethanol with 2 changes for 5 minutes each
  - c. Immerse in 90% ethanol for 5 minutes
  - d. Immerse in 70% ethanol for 5 minutes
  - e. Immerse in 50% ethanol for 5 minutes
  - f. Immerse in distilled water for 5 minutes
3. Antigen Retrieval (Microwave Method):
  - a. Immerse the slides in a microwave compatible tray containing 10 mM Sodium Citrate buffer (pH 6.0) with 0.05% Tween 20.
  - b. Boil the slides and maintain the sub-boiling temperature for 5 minutes in the microwave. Thereafter, take out the tray very carefully and cool it at room temperature (RT) for about 30 minutes.
  - c. Wash the slides 3 times, 3 minutes each by immersing them in TBST (Tris Buffered Saline having 0.05% Tween 20).
4. Quenching of Endogenous Peroxidase:
  - a. Incubate the slides in 3% hydrogen peroxide prepared in methanol for 15 minutes (at RT, in dark conditions).
  - b. Wash the slides in TBST 3 times, 3 minutes each.
5. Protein Blocking:
  - a. Incubate the sections with background sniper solution at RT for 15 minutes (Biocare Medicals, USA).
  - b. Wash the sections 3 times, 3 min each by immersing the slides in TBST.
6. Primary Antibody:
  - a. Dilute the primary antibody at 5ug/ml concentration using PBS as a diluent.
  - b. Incubate the sections with diluted primary antibody for 90 minutes at RT in a humidified chamber.
  - c. Thereafter, wash the slides 4 times, 5 minutes each with TBST.
7. Probe (Secondary Reagent):
  - a. Incubate with MACH 1 Mouse probe for 15 minutes at RT.
  - b. Incubate for 30 min at room temperature with HRP-Polymer (Biocare Medical, USA).
  - c. Wash the slides with TBST 4 times, 5 minutes each
8. Chromogen:
  - a. Mix 32ul of DAB Chromogen with 1 ml of DAB substrate buffer (Biocare Medical, USA).
  - b. Apply 200ul DAB mixture/section and incubate at RT in dark conditions (few seconds - 5 minutes).
  - c. As soon as an appropriate color develops, rinse the slides with deionized water (2-3 brief rinses).
9. Counter stain:
  - a. Counter stain with Hematoxylin for 30 seconds (Vector Labs, USA).
  - b. Wash in deionized water for 1-2 minutes to clear the extra stain.
  - c. Incubate the slides in bluing solution or Scott's water twice for 2 minutes each time.
10. Dehydrate the sections in increasing grades of alcohols:
  - a. 50% alcohol for 1 minute
  - b. 70% for 1 minute
  - c. 90% for 1 minute
  - d. 95% for 1 minute
  - e. 100% for 1 minute
  - f. Xylene with 2 changes for 2 minutes each
11. Mount with DPX mount and cover-slip glass (Fisher Scientific, USA), carefully not allowing any air bubbles to enter.

**NOTE:-** This protocol is provided as a reference tool only. Depending upon the type of tissues /tissue processing and reagents employed, the end user will need to optimize the final conditions for achieving an expected staining.



### **Novus Biologicals USA**

10730 E. Briarwood Avenue  
Centennial, CO 80112  
USA  
Phone: 303.730.1950  
Toll Free: 1.888.506.6887  
Fax: 303.730.1966  
nb-customerservice@bio-techne.com

### **Bio-Techne Canada**

21 Canmotor Ave  
Toronto, ON M8Z 4E6  
Canada  
Phone: 905.827.6400  
Toll Free: 855.668.8722  
Fax: 905.827.6402  
canada.inquires@bio-techne.com

### **Bio-Techne Ltd**

19 Barton Lane  
Abingdon Science Park  
Abingdon, OX14 3NB, United Kingdom  
Phone: (44) (0) 1235 529449  
Free Phone: 0800 37 34 15  
Fax: (44) (0) 1235 533420  
info.EMEA@bio-techne.com

### **General Contact Information**

www.novusbio.com  
Technical Support: nb-technical@bio-techne.com  
Orders: nb-customerservice@bio-techne.com  
General: novus@novusbio.com

### **Products Related to NBP2-45269**

---

NBP2-33376H	Blue Marker Antibody (6F4-F6) [HRP]
HAF007	Goat anti-Mouse IgG Secondary Antibody [HRP]
NB7539	Goat anti-Mouse IgG (H+L) Secondary Antibody [HRP]
NBP1-43317-0.5mg	Mouse IgG2b Kappa Light Chain Isotype Control (MG2b)

---

### **Limitations**

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit [www.novusbio.com/guarantee](http://www.novusbio.com/guarantee)

Earn gift cards/discounts by submitting a review: [www.novusbio.com/reviews/submit/NBP2-45269](http://www.novusbio.com/reviews/submit/NBP2-45269)

Earn gift cards/discounts by submitting a publication using this product:  
[www.novusbio.com/publications](http://www.novusbio.com/publications)

