

Product Datasheet

CD57 Antibody (IHC539)

NBP2-75981-0.1ml

Unit Size: 0.1 ml

Store at 4C in the dark. Do not freeze.

www.novusbio.com



technical@novusbio.com

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:
www.novusbio.com/NBP2-75981

Updated 9/9/2025 v.20.1

Earn rewards for product
reviews and publications.

Submit a publication at www.novusbio.com/publications

Submit a review at www.novusbio.com/reviews/destination/NBP2-75981



NBP2-75981-0.1ml

CD57 Antibody (IHC539)


Product Information	
Unit Size	0.1 ml
Concentration	Please see the vial label for concentration. If unlisted please contact technical services.
Storage	Store at 4C in the dark. Do not freeze.
Clonality	Monoclonal
Clone	IHC539
Preservative	0.1% Sodium Azide
Isotype	IgM Kappa
Purity	Protein L purified
Buffer	Tris Buffer, pH 7.3 - 7.7, with 1% BSA
Target Molecular Weight	110 kDa

Product Description	
Description	Novus Biologicals Mouse CD57 Antibody (IHC539) (NBP2-75981) is a monoclonal antibody validated for use in IHC. All Novus Biologicals antibodies are covered by our 100% guarantee.
Host	Mouse
Gene ID	27087
Gene Symbol	B3GAT1
Species	Human
Immunogen	Recombinant Human CD57

Product Application Details	
Applications	Immunohistochemistry-Paraffin, Immunohistochemistry
Recommended Dilutions	Immunohistochemistry 1:100 - 1:200, Immunohistochemistry-Paraffin 1:100 - 1:200
Application Notes	Visualization of the antigen present in tissue sections is accomplished in a multi-step immunohistochemical staining process, in conjunction with a horseradish peroxidase (HRP) or alkaline phosphatase (AP) linked detection system. Positive control: Tonsil. The 7mL size is a pre-diluted size and no additional dilutions are required before using this item for the intended application.

Images

Immunohistochemistry-Paraffin: CD57 Antibody (IHC539) [NBP2-75981]
- Staining using CD57 antibody (IHC539) in tonsil

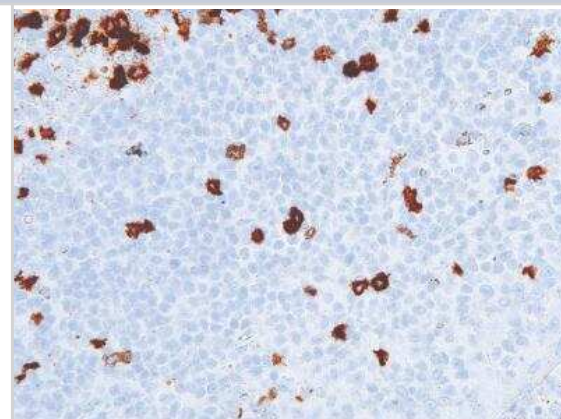

Procedures

Immunohistochemistry Protocol for CD57 Antibody (NBP2-75981)

Immunohistochemistry Protocol

Specimen Collection and Preparation for Analysis

Each tissue section should be fixed with 10% neutral buffered formalin, cut to the applicable thickness (4µm), and placed on a glass slide that is positively charged. The prepared slide may then be baked for a minimum of 30 minutes in a 53-65 degrees C oven (do not exceed 24 hours).

Recommended Staining Protocols

Manual Use:

1. Pretreatment: Perform heat-induced epitope retrieval (HIER) at pH 9 for 10 to 30 minutes.
2. Peroxide Block: Block in peroxidase blocking solution for 5 minutes at room temperature. (Not required if using Alkaline Phosphatase System.)
3. Primary Antibody: Apply antibody directly (Predilute) or dilute antibody at 1:100-1:200 (Concentrate) before applying. Incubate antibody for 10 to 30 minutes at room temperature.
4. Secondary Antibody: Incubate for 20 to 30 minutes at room temperature.
5. Substrate Development: Incubate DAB or Fast Red for 5 to 10 minutes at room temperature.
6. Counterstain: Counterstain with hematoxylin for 0.5 to 5 minutes, depending on the hematoxylin used. Rinse with distilled water and blueing solution for 30 seconds.
7. Dehydrate and apply coverslip.

Automated Staining System:

The stated primary antibody has been optimized and validated using the BOND-MAX fully automated IHC & ISH stainer manufactured by Leica Biosystems, applying IHC Protocol F. The following edits are recommended for the protocol:

- a) Marker Incubation Time: 30 minutes
- b) Heat-induced epitope retrieval (HIER) is recommended using Leica Bond ER Solution 2 for 30 minutes.
- c) Move Peroxide Block step to after Polymer and before Mixed DAB Refine. For all other automated IHC staining systems, refer to the corresponding user manual for specific instructions.





Novus Biologicals USA

10730 E. Briarwood Avenue
Centennial, CO 80112
USA
Phone: 303.730.1950
Toll Free: 1.888.506.6887
Fax: 303.730.1966
nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave
Toronto, ON M8Z 4E6
Canada
Phone: 905.827.6400
Toll Free: 855.668.8722
Fax: 905.827.6402
canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane
Abingdon Science Park
Abingdon, OX14 3NB, United Kingdom
Phone: (44) (0) 1235 529449
Free Phone: 0800 37 34 15
Fax: (44) (0) 1235 533420
info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com
Technical Support: nb-technical@bio-techne.com
Orders: nb-customerservice@bio-techne.com
General: novus@novusbio.com

Products Related to NBP2-75981-0.1ml

HAF007	Goat anti-Mouse IgG Secondary Antibody [HRP]
NB7539	Goat anti-Mouse IgG (H+L) Secondary Antibody [HRP]
NBP1-96975-0.5mg	Mouse IgM Kappa Light Chain Isotype Control (MMK)
210-TA-005	TNF-alpha [Unconjugated]

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NBP2-75981

Earn gift cards/discounts by submitting a publication using this product:
www.novusbio.com/publications

