

# Product Datasheet

## Melanoma Associated Antigen (KBA.62) Antibody (IHC062)

**NBP2-75984-0.1ml**

Unit Size: 0.1 ml

Store at 4C in the dark. Do not freeze.

[www.novusbio.com](http://www.novusbio.com)



[technical@novusbio.com](mailto:technical@novusbio.com)

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:  
[www.novusbio.com/NBP2-75984](http://www.novusbio.com/NBP2-75984)

Updated 9/9/2025 v.20.1

Earn rewards for product  
reviews and publications.

Submit a publication at [www.novusbio.com/publications](http://www.novusbio.com/publications)

Submit a review at [www.novusbio.com/reviews/destination/NBP2-75984](http://www.novusbio.com/reviews/destination/NBP2-75984)



**NBP2-75984-0.1ml**

Melanoma Associated Antigen (KBA.62) Antibody (IHC062)


Product Information	
Unit Size	0.1 ml
Concentration	Please see the vial label for concentration. If unlisted please contact technical services.
Storage	Store at 4C in the dark. Do not freeze.
Clonality	Monoclonal
Clone	IHC062
Preservative	0.1% Sodium Azide
Isotype	IgG1
Purity	Protein A purified
Buffer	Tris Buffer, pH 7.3 - 7.7, with 1% BSA

Product Description	
Description	Novus Biologicals Mouse Melanoma Associated Antigen (KBA.62) Antibody (IHC062) (NBP2-75984) is a monoclonal antibody validated for use in IHC. All Novus Biologicals antibodies are covered by our 100% guarantee.
Host	Mouse
Species	Human
Immunogen	Human KAL Cells derived from Lymphs

Product Application Details	
Applications	Immunohistochemistry-Paraffin, Immunohistochemistry
Recommended Dilutions	Immunohistochemistry 1:100 - 1:200, Immunohistochemistry-Paraffin 1:100 - 1:200
Application Notes	Visualization of the antigen present in tissue sections is accomplished in a multi-step immunohistochemical staining process, in conjunction with a horseradish peroxidase (HRP) or alkaline phosphatase (AP) linked detection system. Positive Control: Melanoma. The 7mL size is a pre-diluted size and no additional dilutions are required before using this item for the intended application.

**Images**

Immunohistochemistry-Paraffin: Melanoma Associated Antigen (KBA.62) Antibody (IHC062) [NBP2-75984] - Staining of Melanoma Associated Antigen (KBA.62) (IHC062) in melanoma.

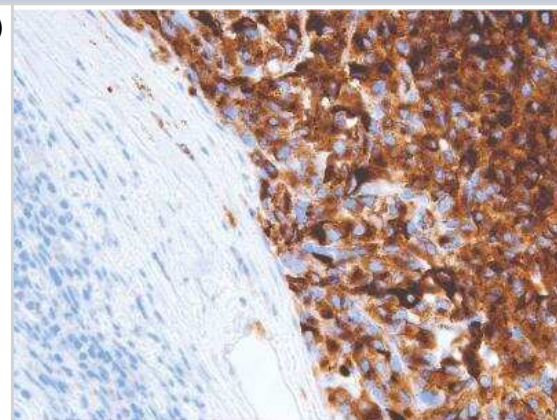

## Procedures

### **Immunohistochemistry protocol for Melanoma Associated Antigen (KBA.62) Antibody (NBP2-75984)**

Immunohistochemistry Protocol for Melanoma Associated Antigen (KBA.62) Antibody (NBP2-75984):

Immunohistochemistry Protocol

#### Specimen Collection and Preparation for Analysis

Each tissue section should be fixed with 10% neutral buffered formalin, cut to the applicable thickness (4um), and placed on a glass slide that is positively charged. The prepared slide may then be baked for a minimum of 30 minutes in a 53-65 degrees C oven (do not exceed 24 hours).

#### Recommended Staining Protocols

##### Manual Use:

1. Pretreatment: Perform heat-induced epitope retrieval (HIER) at pH 9 for 10 to 30 minutes.
2. Peroxide Block: Block in peroxidase blocking solution for 5 minutes at room temperature. (Not required if using Alkaline Phosphatase System.)
3. Primary Antibody: Apply antibody directly (Predilute) or dilute antibody at 1:100-1:200 (Concentrate) before applying. Incubate antibody for 10 to 30 minutes at room temperature.
4. Secondary Antibody: Incubate for 20 to 30 minutes at room temperature.
5. Substrate Development: Incubate DAB or Fast Red for 5 to 10 minutes at room temperature.
6. Counterstain: Counterstain with hematoxylin for 0.5 to 5 minutes, depending on the hematoxylin used. Rinse with distilled water and blueing solution for 30 seconds.
7. Dehydrate and apply coverslip.

##### Automated Staining System:

The stated primary antibody has been optimized and validated using the BOND-MAX fully automated IHC & ISH stainer manufactured by Leica Biosystems, applying IHC Protocol F. The following edits are recommended for the protocol:

- a) Marker Incubation Time: 30 minutes
- b) Heat-induced epitope retrieval (HIER) is recommended using Leica Bond ER Solution 2 for 30 minutes.
- c) Move Peroxide Block step to after Polymer and before Mixed DAB Refine. For all other automated IHC staining systems, refer to the corresponding user manual for specific instructions.





### **Novus Biologicals USA**

10730 E. Briarwood Avenue  
Centennial, CO 80112  
USA  
Phone: 303.730.1950  
Toll Free: 1.888.506.6887  
Fax: 303.730.1966  
nb-customerservice@bio-techne.com

### **Bio-Techne Canada**

21 Canmotor Ave  
Toronto, ON M8Z 4E6  
Canada  
Phone: 905.827.6400  
Toll Free: 855.668.8722  
Fax: 905.827.6402  
canada.inquires@bio-techne.com

### **Bio-Techne Ltd**

19 Barton Lane  
Abingdon Science Park  
Abingdon, OX14 3NB, United Kingdom  
Phone: (44) (0) 1235 529449  
Free Phone: 0800 37 34 15  
Fax: (44) (0) 1235 533420  
info.EMEA@bio-techne.com

### **General Contact Information**

www.novusbio.com  
Technical Support: nb-technical@bio-techne.com  
Orders: nb-customerservice@bio-techne.com  
General: novus@novusbio.com

### **Products Related to NBP2-75984-0.1ml**

---

HAF007	Goat anti-Mouse IgG Secondary Antibody [HRP]
NB7539	Goat anti-Mouse IgG (H+L) Secondary Antibody [HRP]
NBP1-97005-0.5mg	Mouse IgG1 Isotype Control (MG1)

---

### **Limitations**

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit [www.novusbio.com/guarantee](http://www.novusbio.com/guarantee)

Earn gift cards/discounts by submitting a review: [www.novusbio.com/reviews/submit/NBP2-75984](http://www.novusbio.com/reviews/submit/NBP2-75984)

Earn gift cards/discounts by submitting a publication using this product:  
[www.novusbio.com/publications](http://www.novusbio.com/publications)

