

# Product Datasheet

## Exosome Standards (A549 cell line)

### NBP3-11645-200ug

Unit Size: 2 x 100ug Vials

Store at 4C. After reconstitution store at -70C.

[www.novusbio.com](http://www.novusbio.com)



[technical@novusbio.com](mailto:technical@novusbio.com)

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:  
[www.novusbio.com/NBP3-11645](http://www.novusbio.com/NBP3-11645)

Updated 4/13/2025 v.20.1

Earn rewards for product  
reviews and publications.

Submit a publication at [www.novusbio.com/publications](http://www.novusbio.com/publications)

Submit a review at [www.novusbio.com/reviews/destination/NBP3-11645](http://www.novusbio.com/reviews/destination/NBP3-11645)



**NBP3-11645-200ug**

Exosome Standards (A549 cell line)

<b>Product Information</b>	
<b>Unit Size</b>	2 x 100ug Vials
<b>Concentration</b>	Please see the protocols for proper use of this product. If no protocol is available, contact technical services for assistance.
<b>Storage</b>	Store at 4C. After reconstitution store at -70C.
<b>Reconstitution Instructions</b>	Add deionized water, 100 ul for Standard 100 ug and 30 ul for Standard 30 ug, to get a final concentration of 1 mg/mL. Resuspend exosomes pipetting the solution up and down 10-15 times, avoiding bubbles. Vortex the reconstituted standard for 60 seconds.
<b>Buffer</b>	Lyophilized from cell culture media


<b>Product Description</b>	
<b>Description</b>	<p>Highly pure, lyophilized exosome standards with superior stability, optimal for multiple applications including:</p> <ul style="list-style-type: none"> <li>• Assay calibration</li> <li>• Spike-in control for exosome quantification</li> <li>• Protein marker analysis for different techniques such as Western Blot and Flow Cytometry</li> <li>• Extraction and analysis of exosomal RNA and DNA</li> <li>• Standardized positive controls for evaluating immunocapture performance</li> </ul> <p>Quantity per vial of 30 ug size (number of particles in 30 ug: &gt; 1x10<sup>8</sup>). Quantity per vial of 100 ug size (number of particles in 100 ug: &gt; 1x10<sup>10</sup>).</p>
<b>Preparation Method</b>	Isolation involves Tangential flow filtration combined with Size Exclusion Chromatography. Exosomes (small EVs) are quantified and validated for protein content and particle number by Nanoparticle Tracking Analysis as well as for common tetraspanins marker validation. Lyophilization does not alter stability of exosome proteins and nucleic acids.

<b>Product Application Details</b>	
<b>Applications</b>	Simple Western, ELISA, Electron Microscopy, Flow Cytometry, Nucleic Acid Extraction
<b>Recommended Dilutions</b>	Simple Western, Flow Cytometry, ELISA, Electron Microscopy, Nucleic Acid Extraction




## Images


Flow Cytometry: Exosome Standards (A549 cell line) [NBP3-11645] - Phenotyping assays by FACS. Reconstituted Exosomes can be used for profiling biomarkers by FACS analysis. Recommended quantity: 5 ug of reconstituted Exosomes Standards for each test.

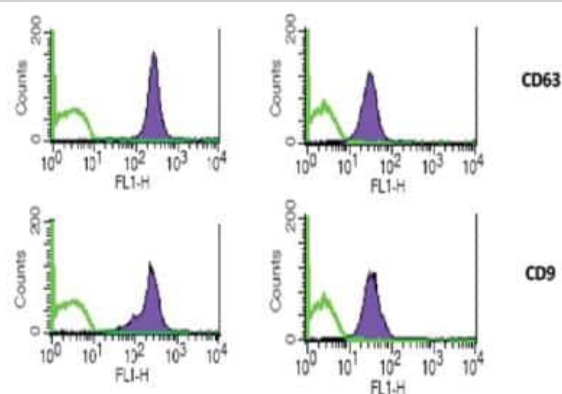

Simple Western: Exosome Standards (A549 cell line) [NBP3-11645] - Simple Western lane view shows lysates of Exosome Standards (A549 cell line), exosome standards (HT29 cell line), and COLO 205 human colorectal adenocarcinoma cell line whole cell lysate (WCL), loaded at 0.2 mg/mL. A specific band was detected for Annexin V at approximately 38 kDa (as indicated) using 50 ug/mL of Mouse Anti-Human Annexin V Monoclonal Antibody. This experiment was conducted under reducing conditions and using the 12-230 kDa separation system.

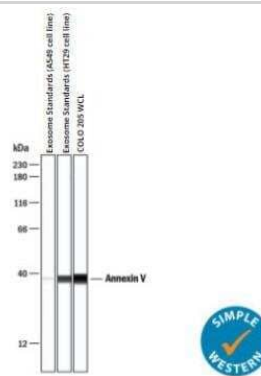

Simple Western: Exosome Standards (A549 cell line) [NBP3-11645] - Simple Western lane view shows lysates of Exosome Standards (A549 cell line), exosome standards (HT29 cell line), and COLO 205 human colorectal adenocarcinoma cell line whole cell lysate (WCL), loaded at 0.2 mg/mL. A specific band was detected for Flotillin-1 at approximately 51 kDa (as indicated) using 1:250 ug/mL of Flotillin-1 Antibody (NBP1-79022). This experiment was conducted under reducing conditions and using the 12-230 kDa separation system.

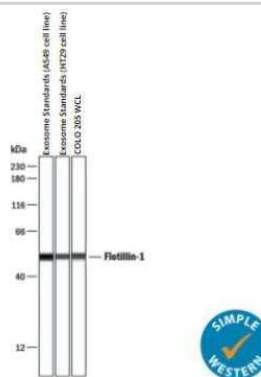



### **Novus Biologicals USA**

10730 E. Briarwood Avenue  
Centennial, CO 80112  
USA  
Phone: 303.730.1950  
Toll Free: 1.888.506.6887  
Fax: 303.730.1966  
nb-customerservice@bio-techne.com

### **Bio-Techne Canada**

21 Canmotor Ave  
Toronto, ON M8Z 4E6  
Canada  
Phone: 905.827.6400  
Toll Free: 855.668.8722  
Fax: 905.827.6402  
canada.inquires@bio-techne.com

### **Bio-Techne Ltd**

19 Barton Lane  
Abingdon Science Park  
Abingdon, OX14 3NB, United Kingdom  
Phone: (44) (0) 1235 529449  
Free Phone: 0800 37 34 15  
Fax: (44) (0) 1235 533420  
info.EMEA@bio-techne.com

### **General Contact Information**

www.novusbio.com  
Technical Support: nb-technical@bio-  
techne.com  
Orders: nb-customerservice@bio-techne.com  
General: novus@novusbio.com

### **Limitations**

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Support products are guaranteed for 6 months from date of receipt.

For more information on our 100% guarantee, please visit [www.novusbio.com/guarantee](http://www.novusbio.com/guarantee)

Earn gift cards/discounts by submitting a review: [www.novusbio.com/reviews/submit/NBP3-11645](http://www.novusbio.com/reviews/submit/NBP3-11645)

Earn gift cards/discounts by submitting a publication using this product:  
[www.novusbio.com/publications](http://www.novusbio.com/publications)

